

FIGURE 1

The diagram illustrates a tissue processing chamber system with the following components and connections:

- Control Device (28):**
 - Processor (54):** Connected to the Interface Board (62) and Memory (58).
 - Interface Board (62):** Connected to the Processor (54).
 - Memory (58):** Contains a **Look-up Table (59)** and **RAM (60)**.
 - Power Board (64):** Contains **Drivers (65)** and is connected to the Processor (54).
- Tissue Processing Chamber (16):**
 - Agitator (36):** For mixing the tissue.
 - Pressure Sensor (38):** Monitors internal pressure.
 - Temperature Sensor (39):** Monitors internal temperature.
 - Fluid Level Sensor (37):** Monitors the fluid level.
 - Valve (34):** Controls fluid flow between the chamber and the paraffin oven.
- Paraffin Oven (32):**
 - Heater (21):** Provides heat for the paraffin.
 - Paraffin Reservoirs (20):** Store the paraffin medium.
- Input/Output Device (30):**
 - Monitor (66):** Displays system status.
 - Keypad (68):** Allows user input.
 - External Ports (70):** For additional connections.

Arrows indicate the flow of control signals and fluid between these components.

1. The first step is to identify the problem. This involves understanding the current situation and the goals that need to be achieved.

Figure 2b


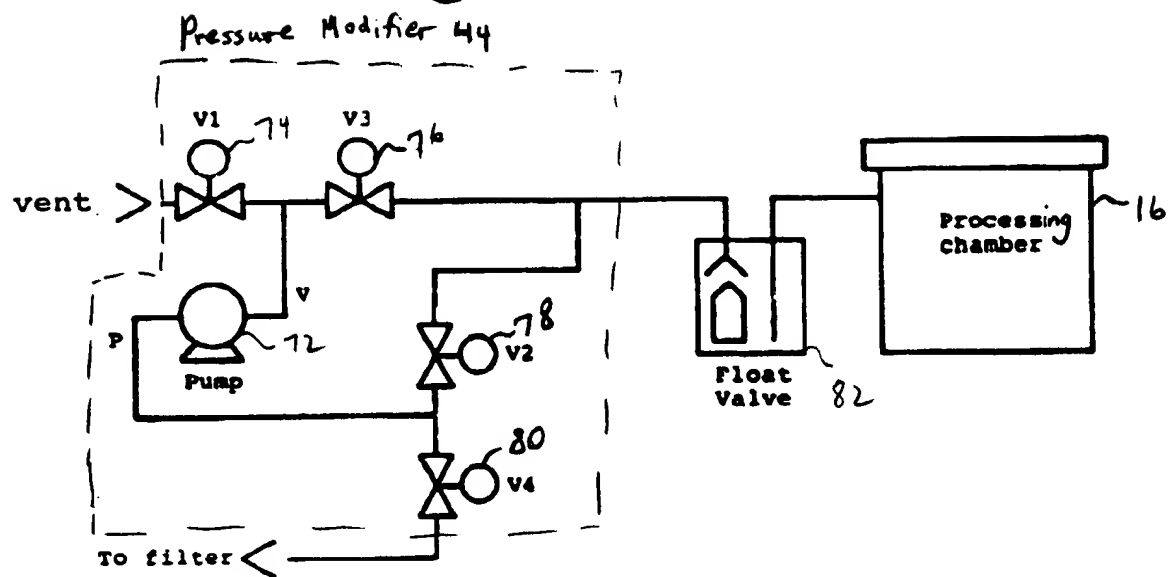


FIGURE 3

2025 RELEASE UNDER E.O. 14176

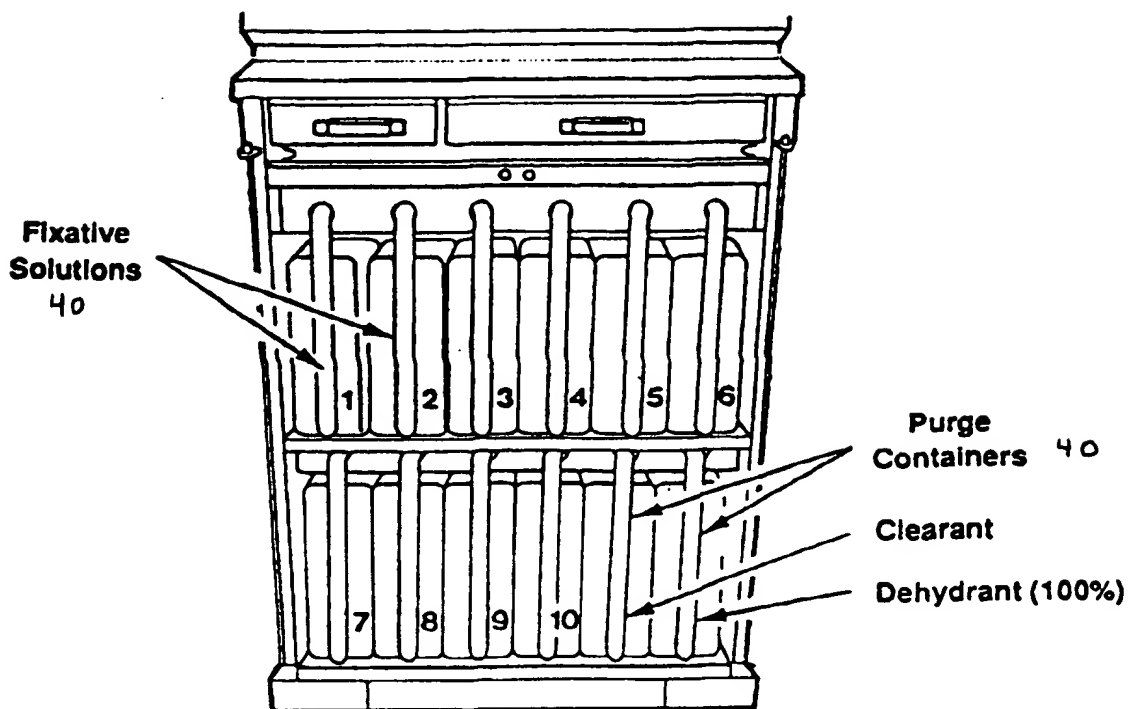

Solvent Container Positions in the Reagent Module

Figure 4


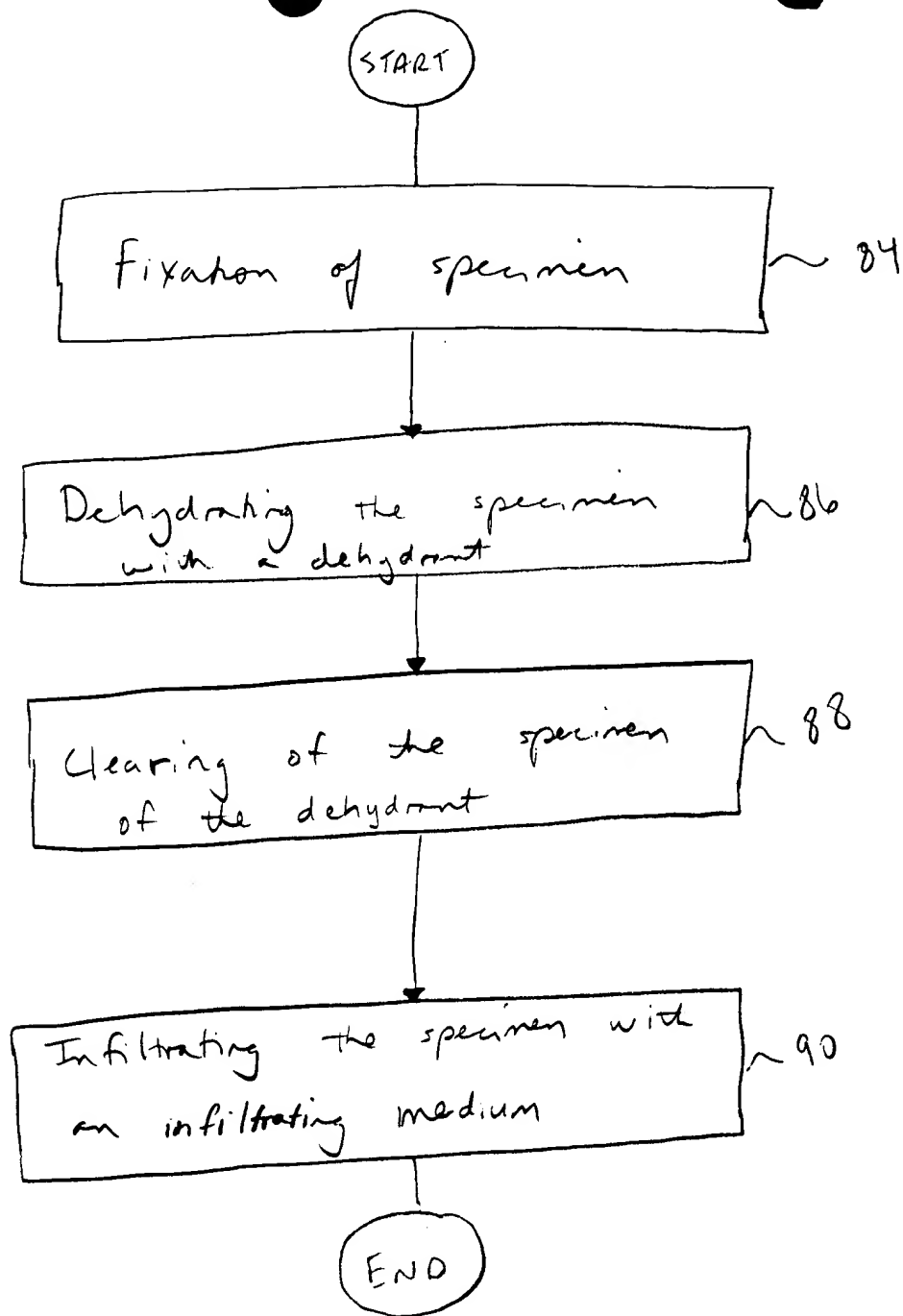


FIGURE 5

FIGURE 6a

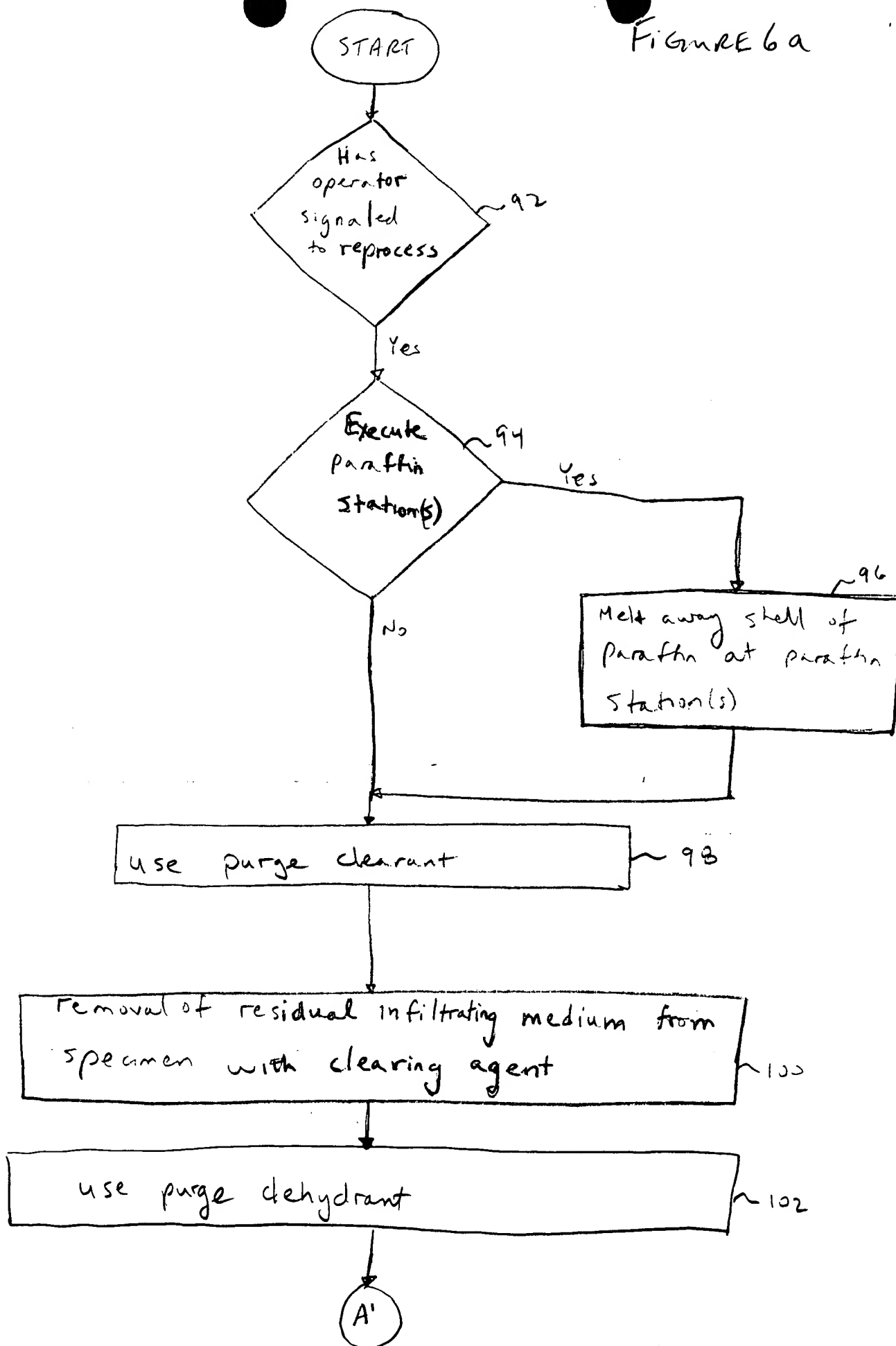

A'

Removal of residual clearing agent by saturating specimen with dehydrating agent ~104

Removal of dehydrating agent with aqueous fluid ~106

Has operator indicated to process specimen ~108

Yes


refixed specimen is dehydrated using dehydrant ~110

dehydrant in specimen is replaced using clearing agent ~112

clearing agent is replaced using infiltrating material ~114

end

FIGURE 6b

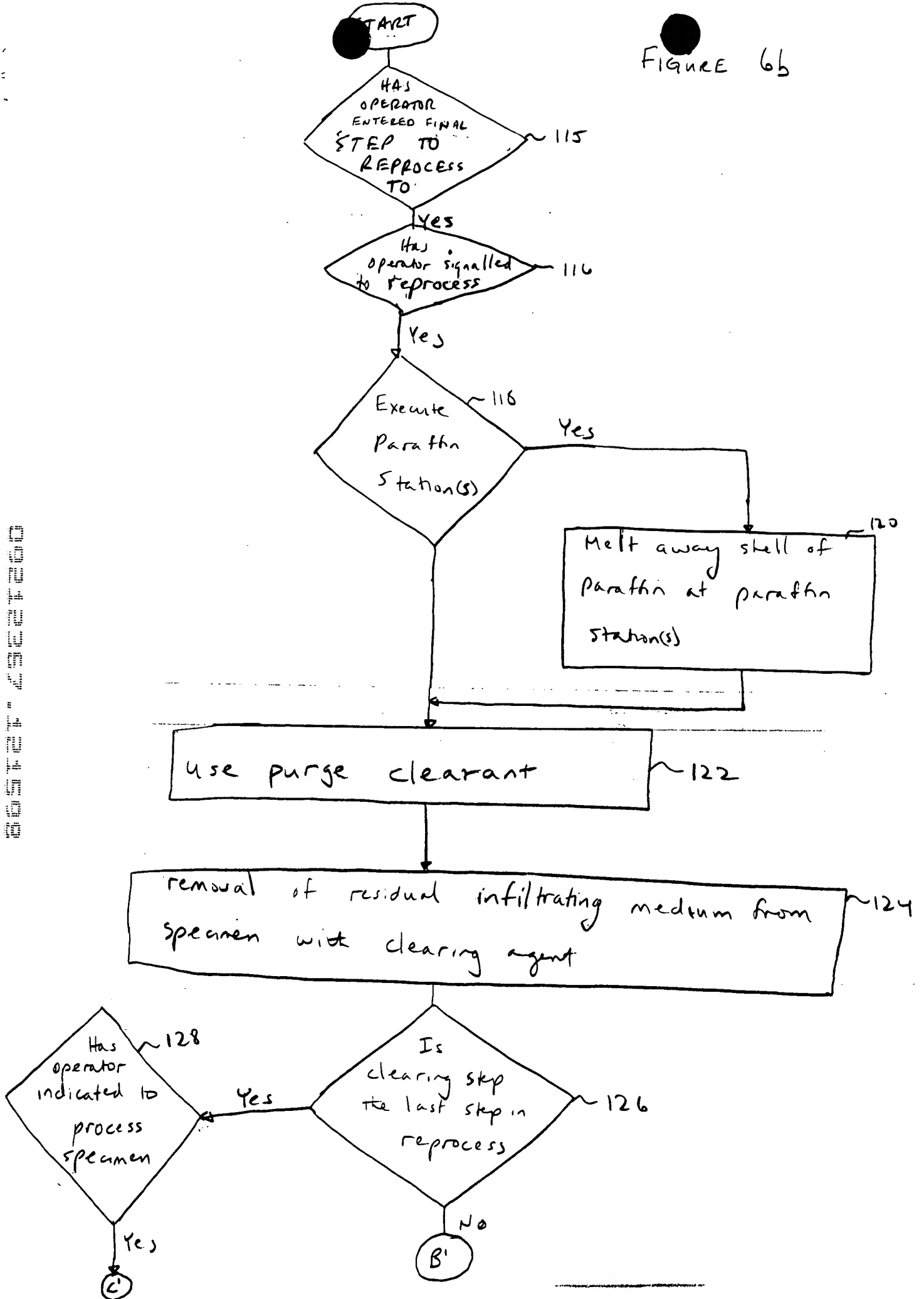

C'

B'

use purge dehydrant ~ 130

removal of residual clearing agent by saturating specimen with dehydrating agent

Has operator indicated to process specimen ~ 136

Is dehydrating step the last step in reprocess ~ 134

Removal of dehydrating agent with aqueous fluid ~ 138

Has operator indicated to process specimen ~ 140

refixed specimen is dehydrated using dehydrant ~ 142

dehydrant in specimen is replaced with clearing agent ~ 144

clearing agent is replaced using infiltrating material ~ 146

END

130
132
134
136
138
140
142
144
146