

UTILITY PATENT APPLICATION TRANSMITTAL

Submit an original and a duplicate for fee processing
(Only for new nonprovisional applications under 37 CFR 1.53(b))

ADDRESS TO:

Assistant Commissioner for Patents
Box Patent Application
Washington, D.C. 20231

Attorney Docket No. 98,714

First Named Inventor David Baunoch, et al.

Express Mail No. EM415807655US

Total Pages 65

APPLICATION ELEMENTS

1. ☒ Transmittal Form with Fee
2. ☒ Specification (including claims and abstract) [Total Pages 50]
3. ☒ Drawings [Total Sheets 10]
4. ☐ Oath or Declaration [Total Pages]
 - a. ☐ Newly executed
 - b. ☐ Copy from prior application
[Note Boxes 5 and 17 below]
 - i. ☐ Deletion of Inventor(s) Signed statement attached deleting inventor(s) named in the prior application
5. ☐ Incorporation by Reference: The entire disclosure of the prior application, from which a copy of the oath or declaration is supplied under Box 4b, is considered as being part of the disclosure of the accompanying application and is hereby incorporated by reference therein.
6. ☐ Microfiche Computer Program
7. ☐ Nucleotide and/or Amino Acid Sequence Submission
 - a. ☐ Computer Readable Copy
 - b. ☐ Paper Copy
 - c. ☐ Statement verifying above copies

ACCOMPANYING APPLICATION PARTS

8. ☐ Assignment Papers
9. ☐ Power of Attorney
10. ☐ English Translation Document (if applicable)
11. ☐ Information Disclosure Statement (IDS)
☐ PTO-1449 Form
☐ Copies of IDS Citations
12. ☐ Preliminary Amendment
13. ☒ Return Receipt Postcard (Should be specifically itemized)
14. ☐ Small Entity Statement(s)
☐ Enclosed
☐ Statement filed in prior application; status still proper and desired
15. ☐ Certified Copy of Priority Document(s)
16. ☐ Other:

17. If a **CONTINUING APPLICATION**, check appropriate box and supply the requisite information:
☐ Continuation ☐ Divisional ☐ Continuation-in-part of prior application Serial No.

APPLICATION FEES


BASIC FEE				\$790.00
CLAIMS	NUMBER FILED	NUMBER EXTRA	RATE	
Total Claims	22 -20=	2	x \$22.00	\$44.00
Independent Claims	6 -3=	3	x \$ 82.00	\$246.00
<input type="checkbox"/> Multiple Dependent Claims(s) if applicable			+\$270.00	\$
Total of above calculations =				\$
Reduction by 50% for filing by small entity =				\$()
<input type="checkbox"/> Assignment fee if applicable			+ \$40.00	\$
TOTAL =				\$1080.00

UTILITY PATENT APPLICATION TRANSMITTAL

Attorney Docket No. 98,714

18. ☐ Please charge my Deposit Account No. 13-2490 in the amount of \$.
19. ☒ A check in the amount of \$1080.00 is enclosed.
20. The Commissioner is hereby authorized to credit overpayments or charge any additional fees of the following types to Deposit Account No. 13-2490:
- a. ☒ Fees required under 37 CFR 1.16.
 - b. ☒ Fees required under 37 CFR 1.17.
 - c. ☒ Fees required under 37 CFR 1.18.
21. ☒ The Commissioner is hereby generally authorized under 37 CFR 1.136(a)(3) to treat any future reply in this or any related application filed pursuant to 37 CFR 1.53 requiring an extension of time as incorporating a request therefor, and the Commissioner is hereby specifically authorized to charge Deposit Account No. 13-2490 for any fee that may be due in connection with such a request for an extension of time.

22. CORRESPONDENCE ADDRESS

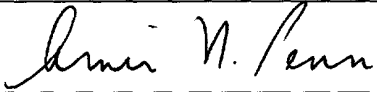Name	McDonnell Boehnen Hulbert & Berghoff
Address	32 nd Floor, 300 South Wacker Drive
City, State, Zip	Chicago, Illinois 60606
23. SIGNATURE OF APPLICANT, ATTORNEY, OR AGENT REQUIRED	
Name	Amir N. Penn, Reg. No. 40,767
Signature	
Date	December 15, 1998

UTILITY (Rev. 11/18/97)

APPLICATION FOR UNITED STATES LETTERS PATENT

SPECIFICATION

(Case No. 98,714)

Title: **METHOD AND APPARATUS FOR AUTOMATED REPROCESSING OF
TISSUE SAMPLES**

Inventors: David Baunoch, a citizen of the United States of America, residing at 11038 N.
Cloud View, Oro Valley, Arizona 85737.

Richard Jones, a citizen of the United States of America, residing at 7110 S.
JenTilly Ln., Tempe, Arizona 85283.

Assignee: Ventana Medical Systems, Inc.

NOTICE REGARDING COPYRIGHT

A portion of the disclosure of this patent document contains matter subject to copyright protection. The copyright owner has no objection to the facsimile reproduction by anyone of the patent disclosure document as it appears in the Patent and Trademark Office files and records but
5 otherwise retains all copyrights whatsoever.

BACKGROUND OF THE INVENTION

A. Field of the Invention

This invention relates generally to the fields of histology and cytology, and more
10 particularly relates to a method and apparatus for reprocessing and processing a specimen.

B. Description of Related Art

Microscopic examination of specimen samples typically involves examining a slice or a cross-section of the sample. In order to obtain a cross-section, the specimen sample undergoes a process to infiltrate the specimen with a paraffin wax or a wax substitute. Thereafter, the block is
15 embedded and sliced into sections using a microtome.

The method of processing the specimen involves fixation of the specimen and preparation of the infiltrated specimen from the fixed specimen. Fixation of the specimen typically involves immersion, subjecting or exposure of the specimen in a fixing agent, such as formalin. Preparation of the infiltrated specimen from the fixed specimen is typically a time-consuming,
20 multi-step process requiring dehydration of the fixed specimen with a dehydrant (such as alcohol), clearing of the dehydrant with a suitable clearant (a typical solvent is xylene), and infiltration of the specimen with an infiltrating medium, such as paraffin wax. In addition, the

dehydration and clearing steps typically require immersion, subjecting or exposure of the specimen in a graded series of reagents for comparatively long periods of time. The time required for tissue preparation may be on the order of 8 to 12 hours. Examples of tissue preparation are in U.S. Patent No. 3,961,097 entitled "Method of Preparing Tissue for Microscopic Examination" and U.S. Patent No. 4,656,047 entitled "Rapid Method for Cell Block Preparation," both of which are hereby incorporated by reference in their entirety.

Different types of specimens, such as any organelle, cell, cell suspension, tissue section, or tissue specimen, may be infiltrated with a paraffin medium for examination. However, different types of specimens may require different types of procedures to be processed properly.

In addition, there may be instances where the specimen may be processed incorrectly, due to contamination of reagents during processing or inadequate fixation. It is typically not until after the specimen has been embedded and sliced that it can be determined whether the specimen has been properly processed. At that point, there are two options: obtain another specimen or reprocess the embedded specimen. If one chooses to reprocess the sample, this involves sequentially immersing, exposing or subjecting the specimen with a series of reagents under controlled conditions. However, this process is very time-consuming and requires a technician to manually proceed through each of the reprocessing steps.

Further, there are instances where a slice or a cross-section of a specimen, after being processed, will be reprocessed for analysis. One instance is ploidy analysis in which tissue sections are cut from the paraffin block, wrapped in a permeable material and reprocessed. The reprocessing steps remove the paraffin using a clearing agent, remove the clearing agent using a dehydrant and remove the dehydrant using an aqueous medium. Nuclei from the specimen are then prepared for DNA analysis using a fluorescent compound.

SUMMARY OF THE INVENTION

In accordance with a first aspect of the invention, an apparatus for automatically reprocessing a specimen from an infiltrated medium to an aqueous fluid is provided. The apparatus has a processing chamber for holding a specimen, means for regulating flow of fluid to the processing chamber, at least one container of a clearant agent, at least one container of a dehydrant agent and at least one container of an aqueous fluid, the containers of clearant, dehydrant and aqueous fluid being connected to the processing chamber via means for regulating flow of fluid to the processing chamber, and a control device having a processor and a memory device, the processor controlling the means for regulating flow of fluid in order to automatically and sequentially connect the processing chamber with the container of clearant agent, the container of dehydrant agent and the container of aqueous solution in order to reprocess the specimen.

In accordance with a second aspect of the invention, a method for automatically reprocessing a specimen using a specimen reprocessing machine having processor for controlling the exposure of the specimen to a clearing agent, a dehydrating agent and an aqueous fluid is provided. The method includes the step of providing the specimen which is infiltrated with an infiltrating medium, indicating to the specimen reprocessing machine that the specimen is to be reprocessed, exposing the specimen to a clearing agent via the processor to remove the infiltrating medium from the specimen, exposing the specimen to a dehydrating agent via the processor to remove the clearing agent, and exposing the specimen to an aqueous fluid via the processor to remove the dehydrating agent from the specimen.

Accordingly, a goal is to process and reprocess specimens for microscopic examination.

These and other objects, features, and advantages of the present invention are discussed or

apparent in the following detailed description.

McDonnell Boeppen Hulbert & Berghoff

BRIEF DESCRIPTION OF THE DRAWINGS

A presently preferred embodiment of the present invention is described herein with reference to the drawings wherein:

FIG. 1 is a front perspective view of the processing and reprocessing system;

FIG. 2a is block diagram of the Operating Module of the processing and reprocessing system;

FIG. 2b is block diagram of the Reagent Module of the processing and reprocessing system;

FIG. 3 is a block diagram of the pressure modifier, float valve and processing chamber in the Operating Module and Reagent Module of Figures 2a and 2b;

FIG. 4 is a front view of the Reagent Module of Figure 1 with the doors removed;

FIG. 5 is a flow chart of the processing of a specimen;

FIG. 6a is a flow chart of the reprocessing of a specimen until introduction of an aqueous fluid in the specimen and processing of the specimen; and

FIG. 6b is a flow chart of the reprocessing of a specimen, until the step as indicated by the operator, and processing of the specimen.

DETAILED DESCRIPTION OF PREFERRED AND ALTERNATIVE EMBODIMENTS OF THE INVENTION

The processing and reprocessing of tissue is accomplished by sequentially putting the specimen (such as any organelle, cell, cell suspension, tissue section, or tissue specimen) to be processed or reprocessed in contact with, or immersed in, a series of reagents under controlled conditions. The reagents may be divided into three types: paraffin, solvents and aqueous

solvents. The conditions that can be controlled while the tissue is in contact with a reagent can be any combination of heat, pressure, vacuum and agitation.

Referring to Figure 1, there is shown one example of a tissue processing and reprocessing system 10. The tissue processing and reprocessing system may consist of two major components: an Operating Module 12 and a Reagent Module 14. The Operating Module and Reagent Module can be placed side-by-side on a benchtop or stacked for a floor mounted configuration, as shown in Figure 1. The specimen is placed in a processing chamber 16, and reagents are sequentially put into the processing chamber 16 from the solvent containers 24, with excess reagents being collected in the reagent drip pan 26 in case of a malfunction. Paraffin is also introduced into the processing chamber 16 with paraffin reservoirs 20. In this apparatus, the control panel 18 indicates the operation of the system 10 and allows for control of the heat, pressure, vacuum and agitation, which affect the processing chamber 16.

An alternative method of putting the specimen in contact with, or immersed in, the reagents is to have each reagent contained in a separate container, and have a mechanical device, such as a robotic arm and controls, to move the specimen from container to container. With this method, the systems for controlling the heat, pressure, vacuum and agitation can, in any combination, be attached to the individual reagent containers or the device for moving the specimen.

Referring to Figure 2a, there is shown a block diagram of the Operating Module 12. The Operating Module 12 houses the processing chamber 16, the control device 28, Input/Output device 30, and the paraffin oven 32 with three reservoirs 20. The processing chamber 16 inside the Operating Module 12 connects through tubing to valve 34, such as a rotary valve (or other means for regulating flow of a fluid) to the paraffin reservoirs, and through tubing to the reagents

in the Reagent Module. In an alternate embodiment, the means for regulating flow of fluid from the paraffin reservoirs to the processing chamber may be performed by any valve, flap, lid, or plug. The processing chamber has an agitator 36, used when the specimen and reagent require stirring. The agitator 36 may be in the form of a rotating stirring device, a recirculating pump, or any other device that causes the reagent to move with respect to the tissue or the tissue to move with respect to the reagent. In addition, the processing chamber has a pressure sensor 38, used to indicate the pressure in the processing chamber 16 to the pressure modifier 44. As described subsequently, the pressure modifier 44 may be accomplished through mechanical means by applying direct mechanical force to the processing chamber through an aneroid, diaphragm, or other mechanical device. The pressure may also be changed by applying pneumatic pressure or vacuum to the processing chamber (*e.g.*, a compressor (air pump) in the system or an external source of vacuum and/or pressure). This may also be accomplished with a mechanical regulator or by cycling the sources of vacuum or pressure on and off.

The Operating Module 12 also includes the paraffin oven 32. Processing and reprocessing of tissue may include the use of an infiltrating medium such as paraffin. The paraffin is stored in a temperature-controlled container in order to keep the paraffin in a liquid state. The temperature of the paraffin reservoirs 20 can be controlled by applying heat directly to the individual containers or by having the paraffin container(s) in a temperature-controlled chamber (such as an oven 32). The oven 32 maintains the paraffin in a liquid state so the system can draw the paraffin into the processing chamber 16, allowing it to penetrate the samples. The processing chamber 16 connects to the paraffin reservoirs through a heated rotary valve 34, which facilitates paraffin selection. At the proper time in the processing and reprocessing program, the rotary valve 34 permits paraffin from the selected reservoir to flow into the

processing chamber, drawn in under vacuum. During the drain cycle, the valve also selects the proper reservoir for the chamber to empty into.

Referring to Figure 2b, there is shown a block diagram of the Reagent Module 14. Processing and reprocessing may require the use of reagents. The Reagent Module 14 contains reagent containers 40 and is connected to the processing chamber 16 in the Operating Module via a solvent/purge line. In one embodiment, the Reagent Module 14 contains twelve reagent containers: ten solvent containers and two purge containers (as shown in Figure 4). The storage temperature of the solvents typically do not need to be controlled and are therefore stored at room temperature. In this arrangement, there is a means for selecting the specific reagent container 40 to move reagents into the processing chamber 16. In one embodiment, the specific reagent container is selected via a set of two valves, one valve 42 (which is set by the processor 54) in the Reagent Module and the second valve 34 (which is set by the processor 54) in the Operating Module. The valve 42 acts as a means for regulating the flow of fluid (which in a preferred embodiment is a liquid and in alternate embodiments may include a liquid, gas or both liquid and gas) between the container 40 and the processing chamber. The valve 42 selects which solvent container connects to the fluid line going to the Operating Module. Thereafter, valve 34 selects which of the paraffin lines or solvent/purge line is connected to the processing chamber. At the proper time in the processing or reprocessing program, the rotary valves 34, 42 (as set by the processor 54) permit only one solvent to flow through the line into the processing chamber, drawn in under vacuum. During the drain cycle, the processor 54 selects the proper setting of valves 34, 42 (as set by the processor 54) for the proper station in order to permit the chamber to empty under pressure. In an alternative embodiment, a valve or other means for regulating the flow of fluid may be mounted on each individual reagent container and connected

to a common manifold which connects to the processing chamber. In another embodiment, other means for regulating the flow of liquid between the containers and the processing chamber include any valve, flap, lid, or plug.

The Reagent Module 14 also has another line to the Operating Module 12 that modifies the pressure in the processing chamber 16. The pressure is modified in the tissue processing and reprocessing system via a pressure modifier 44. As shown in more detail in Figure 3, the pressure modifier 44 serves as a means for introducing and extracting reagents from the processing chamber 16. In one embodiment, this is accomplished by using a pump 72 and a series of valves 74, 76, 78, 80. The pump 72 and valves 74, 76, 78, 80 are contained in the pressure modifier 44 to direct and control the pressure and vacuum. A pressure sensor 38 senses pressure or vacuum in the processing chamber 16. The pump 72 in the Reagent Module cycles on and off as needed to lift fluids into the processing chamber 16 and to drain fluids to their containers 40. In a first state of operation, when valves V1, V2 (74, 78) are closed and valves V3, V4 (76, 80) are open, the pump 72 acts to create a vacuum in the processing chamber 16. In this manner, liquids are drawn into the chamber. Air from the processing chamber 16 is sent to the filter 50, which is described subsequently. In a second state of operation, when valves V3, V4 (76, 80) are closed and valves V1, V2 (74, 78) are open, the pump 72 acts to create a pressure in the processing chamber 16. In this manner, liquids are expelled from the chamber, draining into their respective containers. Air is sent to the processing chamber from a vent. This allows for filling and draining the processing chamber while the specimens remain stationary in the processing chamber. This also permits placing the specimens under vacuum or pressure cycles while immersed in solvents or paraffin to permit thorough infiltration. A pressure differential is created between the storage container and the processing chamber using the pump, to move the

reagents. Alternatively, the force of gravity may be used to move the reagents or paraffin. Further, as shown in Figure 3, there is a float valve 82 which prevents reagents from the processing chamber 16 to enter the pressure modifier 44 in the event that the fluid level in the processing chamber 16 is too high.

5 The Reagent Module 14 has electrical cables for the pump 72, pneumatic valves, 74, 76, 78, 80, rotary valve 42 and blower 48. The Reagent Module 14 features a ventilation system 46 that uses activated charcoal filters to collect solvent fumes before they can escape into the atmosphere. The processing and reprocessing system design reduces the production of fumes. The system handles the fumes from these sources with a built-in ventilating system that filters
10 the air through activated charcoal granules. The system consists of a blower 48 and two filter sections: one for solvent fumes 50, the other for formaldehyde fumes 52. The blower 48 draws air through the Reagent Module 14, up through the filters, and out the back of the unit. This filter system allows operation of the system without the need for a fume hood, external ventilating system, or exhaust fan.

15 The Operating Module 12 further includes a control device 28. The control device, in one embodiment, may be a general purpose computer. This control device 28 automatically controls and sequences the operation of the heaters, motors, pumps and valves, which are controlled via cables. The control device 28 includes, in a preferred embodiment, a processor 54, and in particular, a Hitachi HD-64180 (Z-80) microcontroller. The control device may also
20 include an electro-mechanical timer, an embedded microprocessor circuit, a programmable logic controller, an external computer, or any combination of the above. The control device 28, in one embodiment, contains memory 56 or other computer readable storage medium, including both random access memory (RAM) 58 and read only memory (ROM) 60 in the form of an erasable

programmable read only memory (EPROM). The EPROM contains the system operating program and the text and screen formats for the display. Referring to Appendix A and incorporated herein by reference, there is listed the software having a set of instructions for reprocessing of a specimen. The software is written in Z-80 assembly programming language and is executed on the Hitachi HD-64180 (Z-80) microcontroller.

The control device 28 reads the temperature (via a temperature sensor 39), pressure (via a pressure sensor 38) and the processing chamber fluid level (via a fluid level sensor 37) through the Interface Board 62 and controls the heaters and motors through the Power Board 64. The Power Board 64 contains the drivers 65 for the motors 72, heaters 21, valves 34, 42 and the stirrer 36.

The control device 28 further communicates with the Input/Output device 30 or other user interface. The Input/Output device 30 includes a control panel 18 featuring a monitor 66 such as liquid crystal display (LCD) for displaying menus, instructions and message. The Input/Output device 30 also includes a keypad 68, such as a numeric keypad and an alpha-numeric keyboard or other means for input such as a mouse. The LCD screen assists in programming and operating the system. Through menus, the screen shows status, guides the operator in writing and running reprocessing programs, and serves a variety of maintenance functions. During processing, the monitor 66 shows where specimens are in the cycle, the time in each station, the solution in that station, temperature, and vacuum or pressure. The Input/Output device 30 further includes external ports 70 for connections to external devices such as a printer or a phone line.

The control device 28, in combination with the Input/Output device 30, gives the system its programming flexibility. The operator can program each of the stations (twelve solvent

stations and three paraffin stations) for process time, temperature, vacuum or pressure. The monitor 66 displays all parameters to help the operator while writing the program. A variety of menus give the operator the flexibility of performing a variety of maintenance and service procedures. A special help function gives on-screen assistance at any time without affecting the present status. The operator can tailor processes to match tissue requirements for different solutions and soaking times as well as a combination of heat, pressure and vacuum.

Processing and Reprocessing Cycles

A processing and reprocessing cycle, in one embodiment, consists of filling the processing chamber with a reagent, processing for a programmable amount of time under conditions of controlled temperature, pressure (or vacuum) and agitation. Then draining the reagent back into its storage container. Examples of cycles of the specimen reprocessing and processing system are the fill cycle, the drain cycle and the process cycle.

Fill Cycle

As described previously, four pneumatic valves V1, V2, V3, V4 (74, 76, 78, 80) and the pump 72 perform these cycles, all under computer control.

Before the fill cycle, the system checks that the paraffin oven 32, processing chamber 16 and rotary valve block 34 are up to the programmed temperature. The system then vents the processing chamber 16 and calibrates the pressure sensor 38. The solvent rotary valve 42 moves to the proper position for the selected station and the processing chamber rotary valve 34 moves to the closed position for that station. The system then sets the solenoid valves for vacuum and starts the pump 72. This verifies that the processing chamber 16 and pressure modifier 44 do not leak. If the processing chamber 16 cannot maintain vacuum, the solenoid valves cycle five times to clear any contamination from the valve seats. The system makes a second attempt to establish

vacuum in the chamber. If the processing chamber 16 still cannot maintain vacuum, the system goes to error standby.

If the system successfully established vacuum, then the system vents the processing chamber 16. The processing chamber rotary valve 34 moves to the open position for the selected station. The system sets the solenoid valves 74, 76, 78, 80 for vacuum and cycles the pump 72 on and off to maintain fill vacuum (4 In. Hg for Stations 1-6, 6 in. Hg for Stations 7-10 and purge stations 15 and 16, 2 in. Hg for paraffin stations 11-13, as shown subsequently in Table 1). The system maintains vacuum until the solution triggers the selected level sensor. The processing chamber rotary valve 34 then closes, and the system vents the processing chamber 16.

Drain Cycle

At the beginning of the drain cycle, the system verifies that the paraffin oven 32, processing chamber 16 and valve block 34 are up to temperature. The system then vents the processing chamber 16 and calibrates the pressure sensor 38 by waiting until there is no change in pressure for 1/4 second. The system then stores the pressure sensor reading as the ambient pressure.

The solvent rotary valve 42 moves to the selected station if the station is a solvent station, then the processing chamber rotary valve 34 moves to the closed position for that station. The system next sets the solenoid valves for pressure 74, 76, 78, 80 and starts the pump 72. Similar to the fill cycle, it does this to verify that the processing chamber 16 and the pressure modifier 44 do not leak. If the processing chamber 16 cannot maintain pressure, the solenoid valves will cycle five times to clear any contamination from the valve seats. The system then makes a second attempt to establish pressure in the processing chamber. If the processing chamber still cannot maintain pressure, the system goes to error standby.

The system then releases pressure. The processing chamber rotary valve 34 moves to the open position for the selected station. The system sets the solenoid valves 74, 76, 78, 80 to pressure and the pump 72 starts cycling on and off to maintain drain pressure (1 psi). The system will maintain drain pressure until the processing chamber 16 can no longer hold pressure, indicating that it is empty (the system senses this by the duration of the pump's running cycle). When the system can no longer hold pressure, it vents the processing chamber 16, then waits five seconds for any remaining fluid to collect in the bottom of the processing chamber 16 and its associated plumbing. The system then sets the solenoid valves 74, 76, 78, 80 to pressure, and turns the pump 72 on for two seconds to clear the processing chamber 16 and plumbing of any remaining fluid. The system then vents the processing chamber 16 to release any remaining pressure.

Process Cycle

A programming option allows the specimen processing and reprocessing system to alternate pressure and vacuum while processing or reprocessing tissue to enhance the infiltration of the tissue samples. During programming, the operator sets the values: up to seven pounds per square inches of pressure and a vacuum of up to fifteen inches of mercury. Before the system begins the pressure cycle, it vents the pump while maintaining pressure in the processing chamber. The system does this so that the pneumatic pump starts with no load. After the system starts the pump 72, it sets the solenoid valves 74, 76, 78, 80 to pressure. The cycle runs for 3 minutes at each setting, alternating between vacuum and pressure.

Processing of Specimen

As one example of the specimen processing and reprocessing system, the reagents are arranged in 15 "stations" (3 paraffin stations and 12 solvent stations).

	Station	Reagent	Concentration	Description
	1	Formalin		Fixative
	2	Formalin		Fixative
5	3	Isopropyl Alcohol	70%	Dehydrant
	4	Isopropyl Alcohol	95%	Dehydrant
	5	Isopropyl Alcohol	95%	Dehydrant
	6	Isopropyl Alcohol	100%	Dehydrant
	7	Isopropyl Alcohol	100%	Dehydrant
	8	Alcohol/Xylene	50/50	Dehydrant
10	9	Xylene		Clearant
	10	Xylene		Clearant
	11	Paraffin		
	12	Paraffin		
	13	Paraffin		
15	14	Xylene		Purge Clearant
	15	Isopropyl Alcohol	100%	Purge Dehydrant

Table 1 - Reagent Stations

Depending on the needs in processing the specimen, any number of stations may be present in the machine. In an alternative embodiment in which the specimen is moved from one container of reagents to the next, there may be 15 such containers, as corresponding to the reagents in Table 1, with as a robotic arm and controls to move the specimen from container to container as necessary.

Referring to Figure 5, there is shown a flow chart of a processing of a specimen. The first step involves fixation of the specimen, as shown at block 84. This typically involves immersing the specimen in Formalin, a fixative. In one embodiment of the invention, the specimen is immersed in, subjected to or exposed to a fixing agent at a station, or a multitude of stations, in a processing machine (*see e.g.*, Stations 1 and 2, as shown in Table 1). However, in processing of the sample, the operator may choose to use both stations, only one station or none of the stations (if the specimen has already been immersed in or exposed to a fixative). The specimen is then dehydrated using a dehydrating agent such as alcohol, as shown at block 86. In one embodiment of the invention, the specimen is dehydrated by immersion in, exposure to or being subjected to a

series of alcohol reagents with increasing concentration (*see e.g.*, Stations 3-8, as shown in Table 1). The operator of the processing machine may design a single, or a series, of exposures to alcohol depending on the amount of water contained in the specimen. Thereafter, the specimen is cleared of the dehydrant using a clearing agent, such as xylene (*see e.g.*, Stations 9 and 10, as shown in Table 1), as shown at block 88. Again, depending on processing needs, the specimen may be immersed in, exposed to, or subjected to a single station or both stations. Thereafter, the specimen is infiltrated with an infiltrating medium such as paraffin (*see e.g.*, Stations 11-13, as shown in Table 1), as shown at block 90.

After the specimen has been processed, the machine should be cleaned in order to minimize contaminants of the reagents upon next use the machine. First, a purge clearant, such as xylene, is used in order to clean the paraffin in the processing chamber and the rotary valve on the processing chamber. Second, a purge dehydrant is used to clean any oily residue, or other contaminants, which may be left in the processing chamber. The purge clearant at station 14 is considered to have more impurities of paraffin and other contaminants than, for example, the clearant at station 9. Further, the purge dehydrant at station 15 is considered to have more impurities of oily residue and other contaminants than, for example, the clearant at station 8.

After the tissue has been processed, it is infiltrated with paraffin, embedded in a paraffin block, and sliced into sections using a microtome. At that point, the operator can determine if the specimen has been processed properly. In one instance, the operator may wish to reprocess the remainder of the sample (*i.e.*, the portion of the specimen which has not been sliced up) until the rehydration of the specimen with an aqueous fluid (to a fixing agent, such as formalin, or to water). If that is the case, the operator indicates, via the control panel 18, that the specimen is to

be reprocessed. In addition, if the operator wishes to reprocess a section of the specimen, such as for ploidy analysis, the operator indicates, via the control panel 18, to reprocess the specimen.

Reprocessing of Specimen

Referring to Figure 6a, there is shown a flow chart of the automatic reprocessing of the tissue until rehydration of the specimen and then processing of the specimen. The system, in one embodiment, may wait until the operator has signaled to reprocess the specimen, as shown at block 92. The software therefore has an initiating routine, waiting until the operator has initiated reprocessing. To reprocess tissue, the infiltrating medium is first removed. Typically, a specimen is not only infiltrated with a medium, but also embedded or encased in the same medium. For example, a specimen may be infiltrated with paraffin, and for ease of slicing, may also be embedded or enveloped with a paraffin shell. To remove the paraffin shell, the operator may simply slice the shell away from the specimen with a knife. Otherwise, the operator may allow the reprocessing machine to remove the shell of paraffin, as shown at block 94. Operator input, via the control panel 18, indicates whether the paraffin station(s) are to be used. This may optionally be done by running processing cycles with one or more paraffin stations. (*see e.g.*, Stations 11-13, as shown in Table 1). The shell of paraffin is removed from the specimen by raising the temperature of the tissue to the melting point of the paraffin that has infiltrated the tissue, as shown at block 96.

If the paraffin stations have already been run, the valve 34 and processing chamber 16 are contaminated with paraffin; therefore, the purge clearant should be used. As described previously, the purge clearant is typically used in cleaning the valve 34 and processing chamber 16 when processing a sample. Thus, the purge clearant is already contaminated with paraffin and may clean the valve 34 and processing chamber 16. As a general matter (even if the paraffin

station(s) have not been run), the order of the clearant stations may optionally be from the most contaminated (with paraffin) to the cleanest. This is due to the fact that in removing the paraffin, the clearing agent may become contaminated. In order to avoid contamination of the "cleaner" clearants, the purge clearant should be used first. Otherwise, the cleaner clearants (such as Stations 9 or 10) would be contaminated with paraffin if used directly after a paraffin step. If that were the case, upon processing of a sample again, the clearant in station 9 or 10 would have to be replaced due to contamination. Therefore, the specimen is immersed in, subjected to or exposed to a purge clearant first, as shown at block 98. (*see e.g.*, Station 14, as shown in Table 1). Typically, the specimen is immersed or exposed to the purge clearant for about 20 minutes with the agitator 36 mixing.

The next step is the removal of residual infiltrating medium from a specimen with "cleaner" clearing agent(s) (an agent that is miscible with the embedding and dehydrating agent), as shown at block 100. This is done by running processing cycles with one or more clearant stations, depending on the needs of reprocessing. (*See e.g.*, Stations 9 and 10, as shown in Table 1). The clearant typically used to remove the paraffin is Xylene.

Again, the valve 34 and the processing chamber 16 may be contaminated with an oily residue left by the clearant. Therefore, the specimen may optionally be immersed in, subjected to or exposed to a purge dehydrant before other dehydrants, as shown at block 102. (*see e.g.*, Station 15, as shown in Table 1). As described previously, the purge dehydrant is typically used in cleaning the valve 34 and processing chamber 16 when processing a sample. Thus, the purge dehydrant is already contaminated. Otherwise, the cleaner dehydrants (such as stations 3-8) would be contaminated with the oily residue if used directly after a clearant step. If that were the case, upon processing of a sample again, the dehydrants would have to be replaced due to

contamination. Therefore, the specimen may be immersed or exposed to the purge dehydrant for about 20 minutes with the agitator 36 mixing.

The next step is the removal of the residual clearing agent by saturating the specimen with a dehydrating agent, as shown at block 104. This step is performed whether or not the purge dehydrant is used. This is accomplished by running processing cycles with dehydrants (typically alcohol) with successively higher concentrations of water in which the specimen is immersed in, subjected to or exposed to dehydrant(s). One or many of the dehydrant stations may be used, depending on the needs of reprocessing. (*See e.g.*, Stations 3-8, as shown in Table 1).

The next step is the removal of the dehydrating agent with an aqueous fluid, as shown at block 106 by the specimen immersed in, exposed to or subjected to the aqueous fluid. The aqueous fluid can be used for storage (such as using an aqueous fluid comprised of water) or used to complete the fixation process prior to the repeating of the specimen processing (such as using an aqueous fluid comprised of a fixative such as formalin).

Optionally, the program may wait to determine if the operator has indicated to process specimen, as shown at block 108. The operator may indicate to process specimen at the beginning of the reprocessing sequence, in the middle of reprocessing, or after reprocessing has completed. Alternatively, the program may immediately processes the specimen without operator input.

The fixed specimen is then processed similar to the process steps of Figure 5. In particular, the refixed specimen is dehydrated using a dehydrant (such as alcohol), as shown at block 110. The dehydrant in the specimen is then replaced using a clearing agent (such as

xylene), as shown at block 112. The clearing agent is then replaced using an infiltrating medium (such as paraffin), as shown at block 114.

Referring to Figure 6b, there is shown a flow chart of an alternate embodiment of the automatic reprocessing of the tissue and then processing of the specimen. The operator, after slicing of the specimen, may be able to determine which step in the previous processing sequence was faulty. For example, if the clearant in the processing sequence was contaminated, upon processing, the clearant may have failed to clear all of the dehydrant, thus leaving the specimen with residual dehydrant. Based on this observation, the operator may enter in the control panel the step to which reprocessing should be done, as shown at block 115. This entry may be stored in a look-up table 59 in RAM 58, so that upon reprocessing, the software may determine which step to reprocess to. Alternatively, the entry in the look-up table may be the step in the processing sequence which was faulty. In the example given above, the step would be the clearing step. In this manner, the reprocessing program may read the entry in the look-up table 59, and stop the reprocessing either at the faulty step or the step prior to the faulty step.

Similar to Figure 6a, the system waits until the operator has signaled to reprocess the specimen, as shown at block 116. The operator may allow the reprocessing machine to remove the shell of paraffin, as shown at block 118. The shell of paraffin is removed from the specimen by raising the temperature of the tissue to the melting point of the paraffin that has infiltrated the tissue, as shown at block 120. If the paraffin stations have already been run, the order of the clearant stations is from the most contaminated (with paraffin) to the cleanest. Therefore, the purge clearant is used first, as shown at block 122. (*see e.g.*, Station 14, as shown in Table 1).

The next step is the removal of residual infiltrating medium from the specimen with a clearing agent (an agent that is miscible with the embedding and dehydrating agent), as shown at

block 124. Thereafter, the look-up table 59 is examined to determine whether the clearing step is the last or final step in the reprocess, as shown at block 126. If so, the program then determines if the operator has indicated to process the specimen, as shown at block 128. If so, the program goes to the infiltrating step, as shown at block 146. If the clearing step is not the last or final step in the reprocess, the purge dehydrant is used before other dehydrants, as shown at block 130. (see e.g., Station 15, as shown in Table 1).

The next step is the removal of the residual clearing agent by saturating or exposing the specimen with a dehydrating agent, as shown at block 132. This step is performed whether or not the purge dehydrant is used. One or many of the dehydrant stations may be used, depending on the needs of reprocessing. (See e.g., Stations 3-8, as shown in Table 1).

Thereafter, the look-up table 59 is examined to determine whether the dehydrating step is the last or final step in the reprocess, as shown at block 134. If so, the program then determines if the operator has indicated to process the specimen, as shown at block 136. If so, the program goes to the clearing step, as shown at block 144. If the dehydrating step is not the last step in the reprocess, the next step is the removal of the dehydrating agent with an aqueous fluid, as shown at block 138. Optionally, the program then waits to determine if the operator has indicated to process specimen, as shown at block 140; otherwise, processing begins without operator input.

The operator may indicate to process specimen at the beginning of the reprocessing sequence, in the middle of reprocessing, or after reprocessing has completed. The fixed specimen is then processed similar to the process steps of Figure 5. In particular, the refixed specimen is dehydrated using a dehydrant (such as alcohol), as shown at block 142. The dehydrant in the specimen is then replaced using a clearing agent (such as xylene), as shown at block 144. The

clearing agent is then replaced using an infiltrating medium (such as paraffin), as shown at block 146.

From the foregoing detailed description, it will be appreciated that numerous changes and modifications can be made to the hardware and software aspects of the invention without departure from the true spirit and scope of the invention. For example, the present invention is not dependent on any specific type of computer architecture or type of protocol. This true spirit and scope of the invention is defined by the appended claims, to be interpreted in light of the foregoing specification.

McDonnell Boeppen Hulbert & Berghoff

Appendix A

*

* REVERSE RUN

*

REVRUN: CALL DISREV

REVR2: CALL GETKEY
CP KBF1 ;START

JP Z,RUNREV

CP KBF2 ;EDIT

JP Z,EDITREV

CP KBF3 ;REPROCESSING PROGRAM

JP Z,REPCRG

CP KBF4 ;CANCEL

JP Z,MAINTN

JR REVR2

RET

*

* DISPLAY PROGRAM - REVERSE

*

DISREV: LD A,(REVSTEP)

PUSH AF

LD A,(STANUB)

PUSH AF

LD A,81

LD (SCRNUB),A

CALL DISPSR

LD A,2

LD (DISROW),A

LD A,3

LD (DISCOL),A

CALL PTSET ;POINT AT PROGRAM

XOR A

INC A

PUSH AF

LD (REVSTEP),A

CALL SOREV

CALL GETRRS

CALL BLKDAT ;BLANK DSPDAT

LD HL,DSPDAT ;BUFFER AREA

LD (HL),00

INC HL

LD A,(STANUB)

McDonnell Boehnen Hulbert & Berghoff

300 South Wacker Drive

Chicago, IL 60606

(312) 913-0001

```

CALL BINASC
INC HL
CALL LINEDT          ;LINE AT A TIME

5    LD  A,(DISROW)
INC  A
LD  (DISROW),A
LD  A,2
LD  (DISCOL),A
10   CALL PTSET          ;POINT AT PROGRAM

LD  HL,DSPDAT
CALL MESOUT

15   POP  AF
CP   15
JR   NZ,DISRE1

20   XOR  A
LD  (EDVAR),A          ;ENABLE PROCESS TIME ROUTINE
CALL REVTP            ;DISPLAY TOTAL PROCESS TIME
CALL DISSS            ;DISPLAY STIR SPEED
CALL DISLS            ;DISPLAY LEVEL SENSOR
CALL DSPRPN           ;DISPLAY REPROCESSING NAME
25
POP  AF
LD  (STANUB),A
POP  AF
LD  (REVSTEP),A
30   RET

*****
*
*   REVERSE RUN - RUN
*
35   *****

RUNREV:  XOR  A
LD  (BEEPFLG),A        ;TURN BEEPER OFF
40
IN0  A,(PORT1C)
RES  5,A                ;DISABLE ALARM RELAY
OUT0 (PORT1C),A

45   CALL ALLOFF
CALL ERRSAV

LD  A,(FLAG3)          ;CHECK IF POT FULL
BIT  6,A                ;
50   JP   NZ,CUSERV      ;IF SO GO TO USER SERVICE

LD  A,84                ;SCREEN 1.2.4
LD  (SCRNUB),A          ;POWER FAIL REENTRY
CALL FUNCT

55   LD  A,1

```

	LD	(REVSTEP),A	
	CALL	ERRCLR	;CLEAR ALL ERRORS.
5	LD	HL,FPCDAT	;SET HOLD TIME TO CURRENT TIME
	LD	DE,HDATE	
	LD	A,4	
	CALL	COPYN	
10	LD	HL,FPCTIM	
	LD	DE,HTIME	
	LD	A,4	
	CALL	COPYN	
15	CALL	ACTREV	;GET FIRST ACTIVE STATION.
	OR	A	;IS ANY ACTIVE?
	JP	NZ,RRUN2	;EXIT IF NO STATIONS ACTIVE
20	LD	A,(REPNUB)	
	OR	A	
	JP	NZ,RRUN14	
	JP	MAINTN	
25	RRUN2:CALL	ACTREV	;GET NEXT ACTIVE STATION.
	LD	A,(REVSTEP)	
	LD	(LASTST),A	
30	CALL	STTITL	;INITIALIZE STATS.
	CALL	ITLVRP	;INIT LEVEL STAT
35	CALL	SOREV	
	LD	A,(STANUB)	
	CALL	GETRRS	
	CALL	SETSTA	
	CALL	FINISH	
40	XOR	A	
	LD	(DSPRMS),A	
	CALL	SPTIME	
45	RRUN3:LD	A,1	;RUN STATE FILL
	LD	(BSTATE),A	
	CALL	DISREV	
50	LD	A,(REVSTEP)	
	CALL	PGMPTR	
	LD	A,(FLAG3)	;THIS IS FOR REENTRY
	BIT	6,A	
	JR	NZ,RRUN6	;IF POT FULL DO NOT RESTART STATION
55	LD	A,(FLAG1)	
	SET	0,A	;INFILTRATION FLAG

	LD	(FLAG1),A	
	CALL	LIDCK	
5	CALL	GETRRS	;GET STATION DATA
	CALL	SETSTA	;SET STATION PARAMETERS.
	CALL	FILPOT	;FILL TISSUE POT
10	LD	A,(STIRSP)	
	LD	(STIRS1),A	
	CALL	SPEED	
15	LD	A,(STANUB)	;CLEAR PURGE FLAG IF PURGE STATIONS RUN
	CP	14	
	JR	NZ,RRUN4	
	LD	A,(PRGFLG)	
	RES	0,A	
	LD	(PRGFLG),A	
20	JR	RRUN5	
	RRUN4 CP	15	
	JR	NZ,RRUN5	
	LD	A,(PRGFLG)	
	RES	1,A	
25	LD	(PRGFLG),A	
	RRUN5 LD	A,6	;INIT V/P FLAGS.
	LD	(VPTIM),A	
30	CALL	STATIT	;INITIALIZE POT STAT.
	RRUN6:CALL	LIDCK	
35	LD	A,2	;RUN STATE PROCESSING
	LD	(BSTATE),A	
	CALL	DISREV	
40	XOR	A	
	LD	(DSPRMS),A	;BLANK RUN MESSAGE.
	CALL	ERRSAV	
45	XOR	A	
	LD	(ESTATE),A	;RESET ERROR STATE
50	LD	A,(STANUB)	
	CALL	GETRRS	
	CALL	SETSTA	
	LD	A,(REVSTEP)	
	CALL	PGMPTR	
	CALL	CLRRX	;CLEAR ALL KEYS ENTERED
55	CALL	STOPST	;ENABLE STOP KEY

```

        LD    A,(FLAG4)
        RES   3,A                ;ENABLE VALVES
        LD    (FLAG4),A

5        LD    A,(DSPFLG)
        RES   0,A                ;DISPLAY ALL FIRST TIME
        LD    (DSPFLG),A

*****

10       RRUN7:CALL INPUT          ;TEST FOR STOP

        CALL STIRON
        CALL POTON
15       CALL VPRUN

        CALL DSPRUN
        CALL STATCHK
        CALL LIDCK

20       LD    A,(CLKFLG)
        BIT   0,A                ;TIMED OUT YET?
        JP    Z,RRUN7

25       *****

        CALL ALLOFF
        CALL STOPUS              ;DISABLE STOP KEY

30       LD    A,(STANUB)
        CALL ERRSAV
        CALL SVSTAT

        LD    A,(REVSTEP)
35       INC   A
        LD    (REVSTEP),A

        CALL ACTREV              ;GET NEXT ACTIVE STATION.
        OR    A
40       JR    Z,RRUN10

        RRUN8:LD A,(LASTST)       ;RECALL LAST
        LD    (REVSTEP),A

45       CALL DISREV

        CALL SOREV
        LD    A,(STANUB)
        CALL GETRRS

50       LD    A,(REVSTEP)
        CALL PGMPTR              ;POINT TO STATION.

        LD    A,3
55       LD    (BSTATE),A        ;RUN STATE DRAIN

```



```

JR      RRUN11

RRUN12:  CALL  ALLOFF

5        LD    A,(STANUB)
        CALL  SVSTAT
        CALL  ERRSAV

        CALL  DRAIN

10       CALL  ERRSAV          ;SAVE ERRORS
        XOR   A
        LD    (BSTATE),A      ;DEACTIVATE RUN

15       LD    A,(PRGFLG)      ;CHECK IF PARAFFINS
        AND   03H              ;HAVE RUN.
        LD    (PRGFLG),A

20       LD    A,(FLAG2)
        RES   7,A
        LD    (FLAG2),A      ;CLEAR LID OPEN FLAG

        JP    MAINTN

25  RRUN13: LD    A,(LASTST)    ;DRAIN AND JUMP TO PROCESSING PROGRAM
        LD    (REVSTEP),A      ;RECALL LAST

        CALL  DISREV

30       CALL  SOREV
        LD    A,(STANUB)
        CALL  GETRRS

        LD    A,(REVSTEP)
35       CALL  PGMPTR          ;POINT TO STATION.

        LD    A,6
        LD    (BSTATE),A      ;RUN STATE DRAIN

40       CALL  DRAIN

        RRUN14: CALL  ERRSAV    ;SAVE ERRORS
        XOR   A
        LD    (BSTATE),A      ;DEACTIVATE RUN

45       LD    A,(PRGFLG)      ;CHECK IF PARAFFINS
        AND   03H              ;HAVE RUN.
        LD    (PRGFLG),A

50       LD    A,(FLAG2)
        RES   7,A
        LD    (FLAG2),A      ;CLEAR LID OPEN FLAG

        LD    A,(REPNUB)
55       LD    (PGMNUB),A
        CALL  GETPGM

```



```

        CALL  FINISH                ;COMPUTE IMEDIATE FINISH TIME
        XOR   A
        LD    (STANUB),A
5        JP    RUN1

*****
*
*      GETS NEXT ACTIVE STATION AND PUTS IT IN REVSTEP. 00H IF NOT. - REV
10      *
*****

ACTREV:  LD    A,(REVSTEP)
        PUSH  AF
15
ACTR1: CALL  SOREV
        LD    A,(STANUB)

        CALL  STONCK
20      OR    A                    ;IS STATION OFF?
        JR    NZ,ACTR3

        LD    A,(REVSTEP)
        CP    15                  ;IS NUMBER BEYOND 15?
25      JP    P,ACTR2
        INC   A                    ;ADVANCE STATION.
        LD    (REVSTEP),A
        JR    ACTR1

30      ACTR2: POP  AF
        LD    (REVSTEP),A
        XOR   A
        RET

35      ACTR3: POP  AF
        RET

*****
*
*      OUTPUT TOTAL PROCESS TIME - REV
40      *
*****

REVTPT:  PUSH  AF
45      PUSH  DE
        PUSH  HL

        LD    A,(EDVAR)
        OR    A
50      JR    NZ,REVT2

        LD    HL,FDATE
        LD    DE,SDATE
        LD    A,6
55      CALL  COPYN

```

```

        LD     DE,CTIM1           ;ZERO ACCUM
        LD     HL,ZERO
        CALL  KOPY

5       XOR    A
        LD     (STAON1),A         ;INIT STATION BITS.
        LD     (STAON2),A

10      LD     A,(STANUB)
        CALL  SAVRRS             ;SAVE STADAT
        PUSH  AF

        LD     A,15
        LD     (STANUB),A

15      REVT1: CALL GETRRS        ;GET NEXT STATION
        LD     BC,TIME1
        LD     HL,TMPR1
        CALL  CFLPT             ;CONVERT TO FLOATING

20      CALL  STAON              ;CHECK IF STATION ON.
        LD     IX,CTIM1
        LD     BC,TMPR1
        LD     HL,CTIM1
        CALL  TIMADD            ;ACCUM TIME

25      LD     A,(STANUB)
        DEC    A
        LD     (STANUB),A
        JR     NZ,REVT1        ;TEST IF 15 STATION

        LD     HL,PTIME
        LD     BC,CTIM1
        LD     A,20H
        CALL  CASCI             ;SAVE TOTAL PROCESS TIME

35      LD     A,18
        LD     (DISROW),A
        LD     A,10
        LD     (DISCOL),A
        CALL  PTSET
        CALL  BLKDAT

40      LD     HL,DSPDAT
        LD     DE,PTIME
        CALL  SFTIME            ;OUTPUT PROCESS TIME
        LD     (HL),1AH

        LD     HL,DSPDAT
        CALL  MESOUT

50      POP    AF
        LD     (STANUB),A
        CALL  GETRRS            ;RESTORE STADAT

55      REVT2 POP    HL

```

```

POP    DE
POP    AF
RET

```

```

5  *****
   *
   *    REVERSE RUN - EDIT
   *
   *****
10  EDITREV: LD  A,82
      LD  (SCRNUB),A
      CALL FUNCT          ;DISPLAY SCREEN

15      LD  A,(FLAG1)
      RES  4,A
      LD  (FLAG1),A      ;CLEAR EDITOR RANGE ERROR

20      LD  A,1
      LD  (REVSTEP),A
      XOR  A
      LD  (EDVAR),A

25  RRE1: LD  A,(REVSTEP)
      CALL SOREV
      CALL GETRRS

30      LD  A,(FLAG1)
      SET  7,A
      LD  (FLAG1),A      ;SET REV VIDEO FLAG

      LD  A,(REVSTEP)
      CALL EDLINE

35  RRE2: CALL GETKEY
      CP   KBUP           ;REV FIELD
      JP   Z,RRFREVP
      CP   KBDN           ;ADV FIELD
      JP   Z,RRFFWD
40      CP   KBF1          ;F1 DONE
      JP   Z,REVRUN
      CP   KBF2           ;F2 SOLUTION EDIT
      JP   Z,REVSOL
      CP   KBF3           ;F3 CHANGE STIR SPEED
45      JP   Z,REVSPD
      CP   KBF4           ;F4 CHANGE LEVEL SENSOR
      JP   Z,REVLEV
      CP   0DH            ;ENTER
      JP   Z,RRDAT
50      CP   KBCL          ;CLEAR ON/OFF
      JR   Z,RONOFF
      AND  07FH
      SUB  030H
      CP   0AH            ;NUMBER
55      JP   P,RRE2
      JP   REVNUB

```

```

RONOFF      CALL  ONOFF
            CALL  GETBAC
            CALL  SAVRRS
5            LD    A,(REVSTEP)
            CALL  EDLINE
            JR     RRE2

```

```

*****
10  *
    *      CHANGE STIR SPEED
    *
*****

```

```

15  REVSPD:   LD     A,(STIRSP)
            INC     A
            CP      11
            JP      M,REVSPD1
            LD      A,0

```

```

20  REVSPD1   LD     (STIRSP),A
            CALL    DISSS
            JP      RRE2

```

```

25  *****
    *
    *      CHANGE LEVEL SENSOR
    *
*****

```

```

30  REVLEV:   LD     A,(LEVFLG)
            XOR     002H           ;TOGGLE
            LD      (LEVFLG),A
            CALL    DISLS
35  JP        RRE2

```

```

*****
    *
    *      CHANGE REPROCESSING PROGRAM
    *
40  *****

```

```

REPCHG:      LD     A,(REPNUB)
            INC     A
45  CP        12
            JP      M,REPCHG1
            LD      A,0

```

```

50  REPCHG1   LD     (REPNUB),A
            CALL    DSPRPN
            JP      REVR2

```

```

*****
    *
55  *      INPUT NUMBER INTO VARIABLE POINTED
    *      TO BY REVSTEP AND EDVAR

```

*

```

5  REVNUB:  PUSH  AF
    CALL  PUTNUB

    LD  A,20H
    CALL FILVAR          ;FILL VARIABLE WITH BLANKS

10  POP  AF
    ADD  A,30H          ;CONVERT TO ASCII.
    PUSH HL
    POP  IY
    LD  (IY+5),A        ;LOAD FIRST CHAR.

15  LD  A,C
    CP  5
    CALL Z,EDDP          ;CHECK IF TIME VARIABLE.
    CP  1
    CALL Z,EDDP          ;CHECK IF TIME VARIABLE.

20  LD  A,(REVSTEP)
    CALL EDLINE          ;OUTPUT VARIABLE.

25  REDN1:CALL GETKEY
    CP  KBF1            ;CANCEL
    JP  Z,RRE1
    CP  0DH             ;ENTER
    JP  Z,RRDAT
    CP  KBCL            ;CLEAR
    JR  Z,AONOFF
    AND  07FH
    SUB  030H
    CP  0AH             ;NUMBER
    JP  P,REDN1
    DEC  B
    JR  Z,REDN2
    ADD  A,30H
    CALL SHFTIN          ;SHIFT IN NEW NUMBER.

40  LD  A,C
    CP  5
    CALL Z,EDDP          ;CHECK IF TIME VARIABLE.
    CP  1
    CALL Z,EDDP          ;CHECK IF TIME VARIABLE.

45  LD  A,(REVSTEP)
    CALL EDLINE          ;OUTPUT VARIABLE.
    JR  REDN1

50  REDN2:INC  B
    JR  REDN1

    AONOFF  CALL  ONOFF
    CALL  GETBAC
55  CALL  SAVRRS
    LD  A,(REVSTEP)

```

CALL EDLINE
JR REDN1

5 *
* SAVES THE FIELD
*

10 RRDAT: CALL RNGCHK
CP 0AAH ;TEST IF IN RANGE.
JP NZ,RRD1
CALL GETRRS
LD A,(FLAG1)
15 SET 4,A
LD (FLAG1),A ;SET EDITOR RANGE ERROR
JP RRE1

RRD1 CALL SAVRRS
20 LD A,(REVSTEP)
CALL OUTLIN
LD A,(FLAG1)
RES 4,A
LD (FLAG1),A ;CLEAR EDITOR RANGE ERROR
25 JP RRFWD

30 *
* REV FIELD
*

RRFREV: LD A,(FLAG1)
RES 7,A
35 LD (FLAG1),A ;CLEAR REV VIDEO FLAG
BIT 4,A ;RETURN IF EDITOR RANGE ERROR
JP NZ,RRE1

40 LD A,(REVSTEP)
CALL EDLINE
CALL REVTPPT ;DISPLAY TOTAL PROGRAM TIME

45 LD A,(EDVAR)
OR A ;FIRST VARIABLE IN STA?
JR Z,RRF1
DEC A
LD (EDVAR),A
JP RRE1

50 RRF1: CALL SAVRRS
LD A,(REVSTEP)
CALL OUTLIN
CP 1 ;FIRST STATION?
55 JR Z,RRF2

```

      DEC  A                      ;ACCESS PREVIOUS LINE
      LD   (REVSTEP),A
5  RRF4: CALL OFFRCK
      OR   A
      JR   NZ,RRF3
      XOR  A
      LD   (EDVAR),A
      JP   RRE1

10  RRF2: LD   A,15
      LD   (REVSTEP),A
      JR   RRF4

15  RRF3: LD   A,3                  ;RAP AROUND
      LD   (EDVAR),A
      JP   RRE1

*****
*
20  *   ADVANCE FIELD
*
*****

RRFFWD: LD   A,(FLAG1)
25      RES  7,A
      LD   (FLAG1),A              ;CLEAR REV VIDEO FLAG

      BIT  4,A                    ;RETURN IF EDITOR RANGE ERROR
30      JP   NZ,RRE1

      LD   A,(REVSTEP)
      CALL EDLINE
      CALL REVTPY                ;DISPLAY TOTAL PROGRAM TIME

35  RAF0 LD   A,(EDVAR)
      CP   3                      ;FIRST VARIABLE IN STA?
      JR   Z,RAF1
      CALL OFFRCK
      OR   A
40      JR   Z,RAF1

      LD   A,(EDVAR)
      INC  A
      LD   (EDVAR),A
45      JP   RRE1

      CALL SAVRRS
      LD   A,(REVSTEP)
50      CALL OUTLIN
      CP   15                    ;FIRST STATION?
      JR   Z,RAF2
      INC  A                      ;ACCESS PREVIOUS LINE
      LD   (REVSTEP),A
      JR   RAF3
55  RAF2: LD   A,1
      LD   (REVSTEP),A

```

```

RAF3: XOR    A                ;RAP AROUND
      LD     (EDVAR),A
      CALL  SOREV
5      CALL  GETRRS
      JP     RRE1

```

OFFRCK: PUSH DE

```

10      CALL  SOREV
      CALL  GETRRS
      LD     DE,TIME1
      CALL  CKZERO

15      POP   DE
      RET

```

```

*
20      *      LOADS STATION DATA IN STANUB TO STADAT - REVERSE
*
*****

```

GETRRS: PUSH AF

```

25      PUSH  BC
      PUSH  HL
      PUSH  DE

      LD     A,(STANUB)
30      CP     16
      JP     P,GSTR1
      LD     HL,STARR
      DEC    A
      CALL  GETPT
35      LD     DE,STADAT
      LD     BC,31
      LDIR

      LD     HL,STADAT
      LD     DE,STADAT1
40      LD     BC,31
      LDIR

```

```

GSTR1: POP   DE
      POP   HL
45      POP   BC
      POP   AF
      RET

```

```

50      *
*      SAVES STATION DATA IN STADAT TO STATION IN STANUB - REVERSE
*
*****

```

```

55      SAVRRS: PUSH AF
      PUSH  BC

```



```

        PUSH HL
        PUSH DE

5         LD  A,(STANUB)
        CP   16
        JP   P,SSTR1
        LD   HL,STARR
        DEC  A
        CALL GETPT
10        PUSH HL
        POP  DE
        LD   HL,STADAT
        LD   BC,31
        LDIR
15        LD   HL,STADAT
        LD   DE,STADAT1
        LD   BC,31
        LDIR

20  SSTR1: POP  DE
        POP  HL
        POP  BC
        POP  AF
        RET

25        *****
        *
        *   EDIT SOLUTION LIST - REVERSE
        *
30        *****

        REVSOL: LD   A,83
                LD   (SCRNUB),A
                CALL DISPSR
35
                CALL ORSOL

        REVSL1: LD   A,(FLAG1)
                SET  7,A           ;REV VIDEO
40        LD   (FLAG1),A

                LD   A,(REVSTEP)
                CALL OUTSOL

45  REVSL2: CALL  GETKEY
        CP   030H           ;SOLUTION LIST
        JP   Z,RSOLST
        CP   KBF1           ;DONE
        JP   Z,REVRUN
50        CP   KBCL           ;CLEAR
        JP   Z,CLRSOL
        CP   KBDN           ;FIELD ADV
        JP   Z,RSADV
        CP   KBUP           ;FIELD REV
55        JP   Z,RSREV
        CP   0DH            ;ENTER

```

```

      JP      Z,REVSL4
      AND     07FH
      SUB     030H
      CP      0AH
5       JP      M,ENRSOL
      JR      REVSL2

```

```

RSOLST: CALL  SOLIST
      JR      REVSOL

```

```

10      REVSL4: CALL  SAVRRS
      JR      REVSL1

```

```

15      CLRSOL: LD   A,0
      LD   (SOLUT1),A
      LD   HL,PERCT1
      LD   A,20H
      CALL FILVAR
      CALL SAVRRS
20      JR      REVSL1

```

```

      ENRSOL: CALL  REVIN
      CALL  SAVRRS
      JR      REVSL1

```

```

25      *****
      *
      *      ADVANCE TO NEXT STATION - REVERSE
      *
30      *****

```

```

      RSADV: LD   A,(FLAG1)
      RES     7,A           ;NO REV VIDEO
      LD   (FLAG1),A
35      LD   A,(REVSTEP)
      CALL  OUTSOL

```

```

      INC     A
40      CP     16
      JR      NZ,RSA1
      LD   A,1
      RSA1: LD   (REVSTEP),A
      CALL  SOREV
45      CALL  GETRRS
      JP      REVSL1

```

```

      *****
      *
50      *      REVIEW PREVIOUS STATION - REVERSE
      *
      *****

```

```

55      RSREV: LD   A,(FLAG1)
      RES     7,A           ;NO REV VIDEO
      LD   (FLAG1),A

```

```

LD    A,(REVSTEP)
CALL  OUTSOL

```

```

5      DEC    A
      JR      NZ,RSR1
      LD      A,15
RSR1:  LD      (REVSTEP),A
      CALL    SOREV
10     CALL    GETRRS
      JP      REVSL1

```

```

*****

```

```

*
```

```

15    *      OUTPUT 15 STATIONS WITH SOLUTIONS. - REVERSE

```

```

*
```

```

*****

```

```

20    ORSOL: PUSH HL
      PUSH    BC

```

```

      LD      A,1
      LD      (REVSTEP),A
      CALL    SOREV
25     CALL    GETRRS

```

```

      LD      A,(FLAG1)
      RES     7,A          ;NO REV VIDEO
      LD      (FLAG1),A

```

```

30    RSOT1: LD      A,(REVSTEP)
      CALL    OUTSOL

```

```

      INC     A
35     CP      16
      JR      Z,RSOT2
      LD      (REVSTEP),A
      CALL    SOREV
      CALL    GETRRS
40     JR      RSOT1

```

```

      RSOT2 LD      A,1
      LD      (REVSTEP),A
      CALL    SOREV
45     CALL    GETRRS

```

```

      POP     BC
      POP     HL
      RET

```

```

50    *****

```

```

*
```

```

*      STATION LIST - REVERSE

```

```

*
```

```

55    *****

```

```

STARR:DW   RR1
          DW   RR2
          DW   RR3
          DW   RR4
5          DW   RR5
          DW   RR6
          DW   RR7
          DW   RR8
          DW   RR9
10         DW   RR10
          DW   RR11
          DW   RR12
          DW   RR13
          DW   RR14
15         DW   RR15

```

```

*****
*
*   STATION ORDER LIST - REVERSE
20  *   GET STANUB FROM REVSTEP
*
*****

```

```

SOREV:PUSH  AF
25          PUSH  BC
          PUSH  HL

          LD      HL,RSTALST
          LD      A,(REVSTEP)
30          DEC   A
          LD      C,A
          LD      B,0
          ADD     HL,BC
          LD      A,(HL)
35          LD      (STANUB),A

          POP     HL
          POP     BC
          POP     AF
40          RET

```

```

RSTALST    DB      13
          DB      12
          DB      11
45          DB      14
          DB      10
          DB      9
          DB      15
          DB      8
50          DB      7
          DB      6
          DB      5
          DB      4
          DB      3
55          DB      2
          DB      1

```

CLAIMS

We claim:

1. An apparatus for automatically reprocessing a specimen from an infiltrating medium to an aqueous fluid comprising in combination:

a processing chamber for holding a specimen;

means for regulating flow of fluid to the processing chamber;

at least one container of a clearant agent, at least one container of a dehydrant agent and at least one container of an aqueous fluid, the containers of clearant, dehydrant and aqueous fluid being connected to the processing chamber via means for regulating flow of fluid to the processing chamber; and

a control device having a processor and a memory device, the processor controlling the means for regulating flow of fluid in order to automatically and sequentially connect the processing chamber with the container of clearant agent, the container of dehydrant agent and the container of aqueous solution in order to reprocess the specimen.

2. The apparatus of claim 1 wherein means for regulating flow of fluid includes a rotary valve and wherein the processor selects the containers of clearant, dehydrant or aqueous fluid by setting the rotary valve.

3. The apparatus of claim 1 further comprising:

at least one container of an infiltrating medium being connected to the processing chamber by a second valve and wherein the processor controls the second valve.

4. The apparatus of claim 3 wherein the processor further controls the means for regulating flow of fluid and the second valve in order to automatically and sequentially, after the connection to the container of aqueous fluid, connect the processing chamber with the container of dehydrant agent, the container of clearant and the container of infiltrating medium in order to process the specimen.

5. The apparatus of claim 1 further comprising a container of purge dehydrant being connected to the processing chamber by the means for regulating flow of fluid, the processor controlling the means for regulating flow of fluid in order to automatically and sequentially connect the processing chamber with the container of clearant agent, the container of purge dehydrant, the container of dehydrant agent and the container of aqueous solution in order to reprocess the specimen.

6. The apparatus of claim 5 further comprising a container of purge clearant being connected to the processing chamber by the means for regulating flow of fluid, the processor controlling the means for regulating flow of fluid in order to automatically and sequentially connect the processing chamber with the container of purge clearant, the container of clearant agent, the container of purge dehydrant, the container of dehydrant agent and the container of aqueous solution in order to reprocess the specimen.

7. A computer readable storage medium containing a set of instructions for a general purpose computer having a user interface comprising means for input, an output driver for

connections to at least one valve, the valve being connected to at least one container of a clearant agent, at least one container of a dehydrant agent and at least one container of an aqueous fluid, the set of instructions comprising:

an initiating routine operatively associated with said user interface for permitting a user to initiate reprocessing via the means for input, said means for input being associated with a reprocessing program accessible to said computer;

a run routine for implementing said reprocessing program selected by the user, the reprocessing program controlling the output drive to the valve in order to automatically and sequentially connect the valve to the container of clearant agent, the container of dehydrant agent and the container of aqueous solution for reprocessing of a specimen.

8. The computer readable storage medium of claim 7 wherein the valve is connected to a container of purge dehydrant and wherein

the reprocessing program automatically and sequentially connects the valve to the container of purge dehydrant after connection of the valve to the container of clearant agent.

9. The computer readable storage medium of claim 8 wherein the valve is connected to a container of purge clearant and wherein

the reprocessing program automatically and sequentially connects the valve to the container of purge clearant before connection of the valve to the container of clearant agent.

10. Method for automatically reprocessing a specimen using a specimen reprocessing machine having processor for controlling the exposure of the specimen to a clearing agent, a

dehydrating agent and an aqueous fluid, the method comprising the steps of:

providing the specimen which is infiltrated with an infiltrating medium;

indicating to the specimen reprocessing machine that the specimen is to be reprocessed;

exposing the specimen to a clearing agent via the processor to remove the infiltrating medium from the specimen; thereafter

exposing the specimen to a dehydrating agent via the processor to remove the clearing agent; and thereafter

exposing the specimen to an aqueous fluid via the processor to remove the dehydrating agent from the specimen.

11. The method of claim 10 wherein the processor further controls the exposure of the specimen to an infiltrating medium and further comprising the steps of:

exposing the specimen to a dehydrating agent via the processor after exposing the specimen to an aqueous fluid;

exposing the specimen to a clearing agent to remove the dehydrating agent; and

exposing the specimen to an infiltrating medium to replace the clearing agent.

12. The method of claim 10 wherein the clearing agent is xylene.

13. The method of claim 10 wherein the dehydrating agent is alcohol.

14. The method of claim 10 wherein the aqueous fluid is formalin.

15. A computer-readable storage medium containing a set of instructions for a general purpose computer, said set of instructions implementing the procedure shown in Figure 6a.

16. Method for reprocessing a specimen which is infiltrated with an infiltrating medium using a specimen reprocessing system having a processing chamber, the method comprising the steps of:

subjecting the specimen to at least one exposure to a clearant;

subjecting the tissue sample to at least one exposure to a purge dehydrant after the exposure to the clearant, the purge dehydrant being contaminated with clearant and being used to clean the processing chamber of the clearant;

subjecting the specimen to at least one exposure to a dehydrant after the exposure to the purge dehydrant; and

subjecting the specimen to an aqueous fluid.

17. The method of claim 16 further comprising the steps of:

subjecting the specimen to at least one exposure to paraffin; and thereafter

subjecting the specimen to at least one exposure to a purge clearant after the exposure to the bath of paraffin, the purge clearant being contaminated with clearant and being used to clean the processing chamber of paraffin.

18. The method of claim 17 further comprising the steps of:

subjecting the specimen to the dehydrant after subjecting the specimen to an aqueous fluid;

subjecting the specimen to the clearant to remove the dehydrant; and
subjecting the specimen to an infiltrating medium to replace the clearant.

19. Method for reprocessing a specimen which is infiltrated with an infiltrating medium using a specimen reprocessing system having a processing chamber, the method comprising the steps of:

inputting a final step in reprocessing the specimen;
subjecting the specimen to at least one exposure to a clearant;
determining whether the step of subjecting the specimen to at least one exposure of clearant is the final step in reprocessing the specimen;
subjecting the specimen to at least one exposure to a dehydrant if the step of subjecting the specimen to at least one exposure of clearant is not the final step in reprocessing the specimen;
determining whether the step of subjecting the specimen to at least one exposure of dehydrant is the final step in reprocessing the specimen; and
subjecting the specimen to an aqueous fluid if the step of subjecting the specimen to at least one exposure of dehydrant is not the final step in reprocessing the specimen.

20. The method of claim 19 further comprising the step of
subjecting the specimen to an infiltrating medium to replace the clearant, after the step of determining whether the step of subjecting the specimen to at least one exposure of clearant is the final step, if the step of subjecting the specimen to at least one exposure of clearant is the final step in reprocessing the specimen.

21. The method of claim 19 further comprising the steps of:

subjecting the specimen to the clearant to remove the dehydrant, after the step of determining whether the step of subjecting the specimen to at least one exposure of dehydrant is the final step, if the step of subjecting the specimen to at least one exposure of dehydrant is the final step in reprocessing the specimen; and

subjecting the specimen to an infiltrating medium to replace the clearant.

22. The method of claim 19 further comprising the steps of:

subjecting the specimen to the dehydrant after subjecting the specimen to an aqueous fluid;

subjecting the specimen to the clearant to remove the dehydrant; and

subjecting the specimen to an infiltrating medium to replace the clearant.

ABSTRACT

A method and apparatus of automatically reprocessing a specimen for microscopic examination is disclosed. Processing of a specimen for microscopic examination involves fixation of the specimen and preparation of the embedded specimen from the fixed specimen. There are instances where, once a specimen has been processed, it is necessary to reprocess the specimen due to contamination of reagents during processing or inadequate fixation. The system automatically reprocesses a specimen by removing residual embedding material from the specimen with a clearing agent, removing the clearing agent with a dehydrating agent, and removing the dehydrating agent with an aqueous fluid.

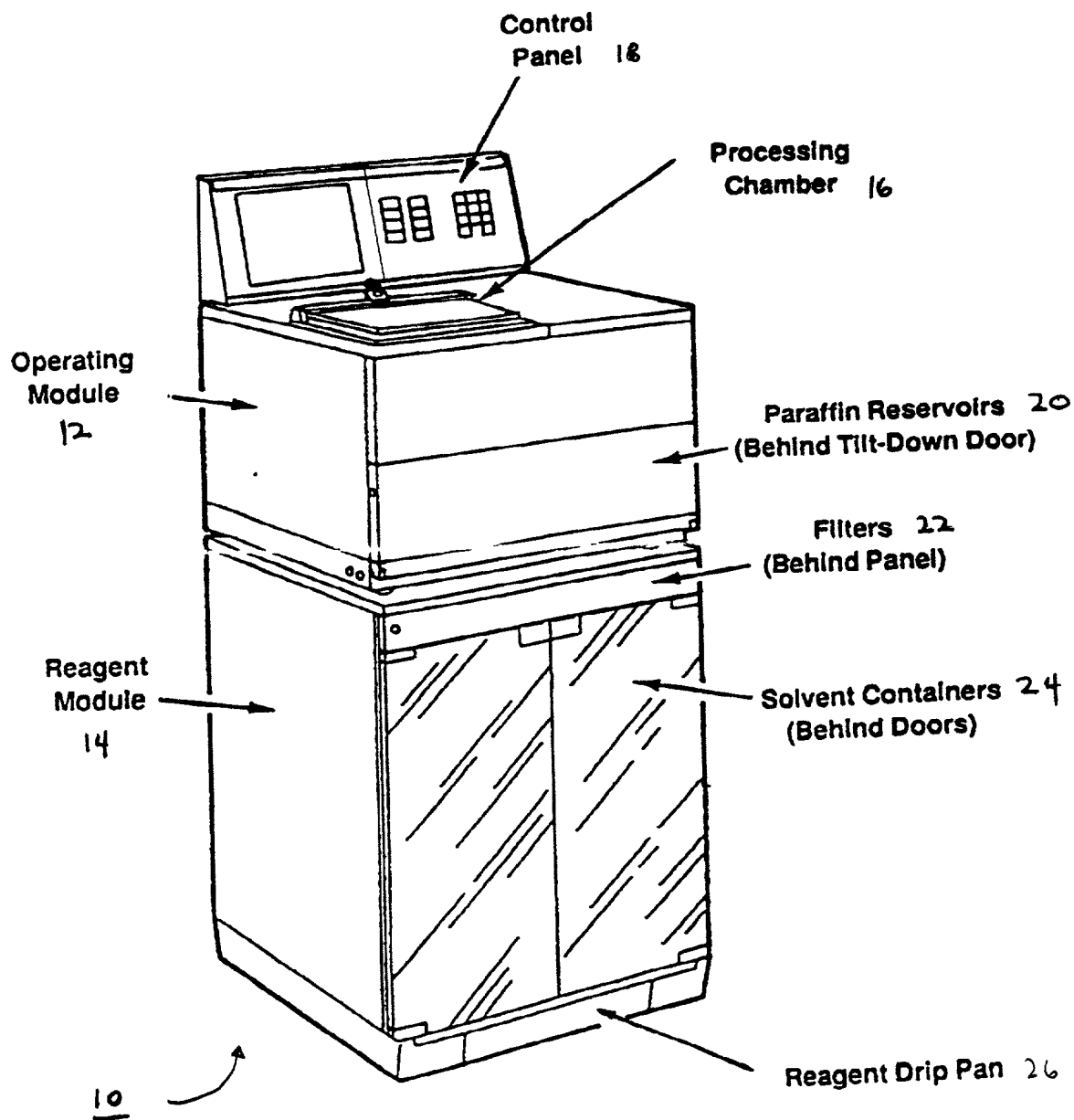


FIGURE 1

Operating Mode 12

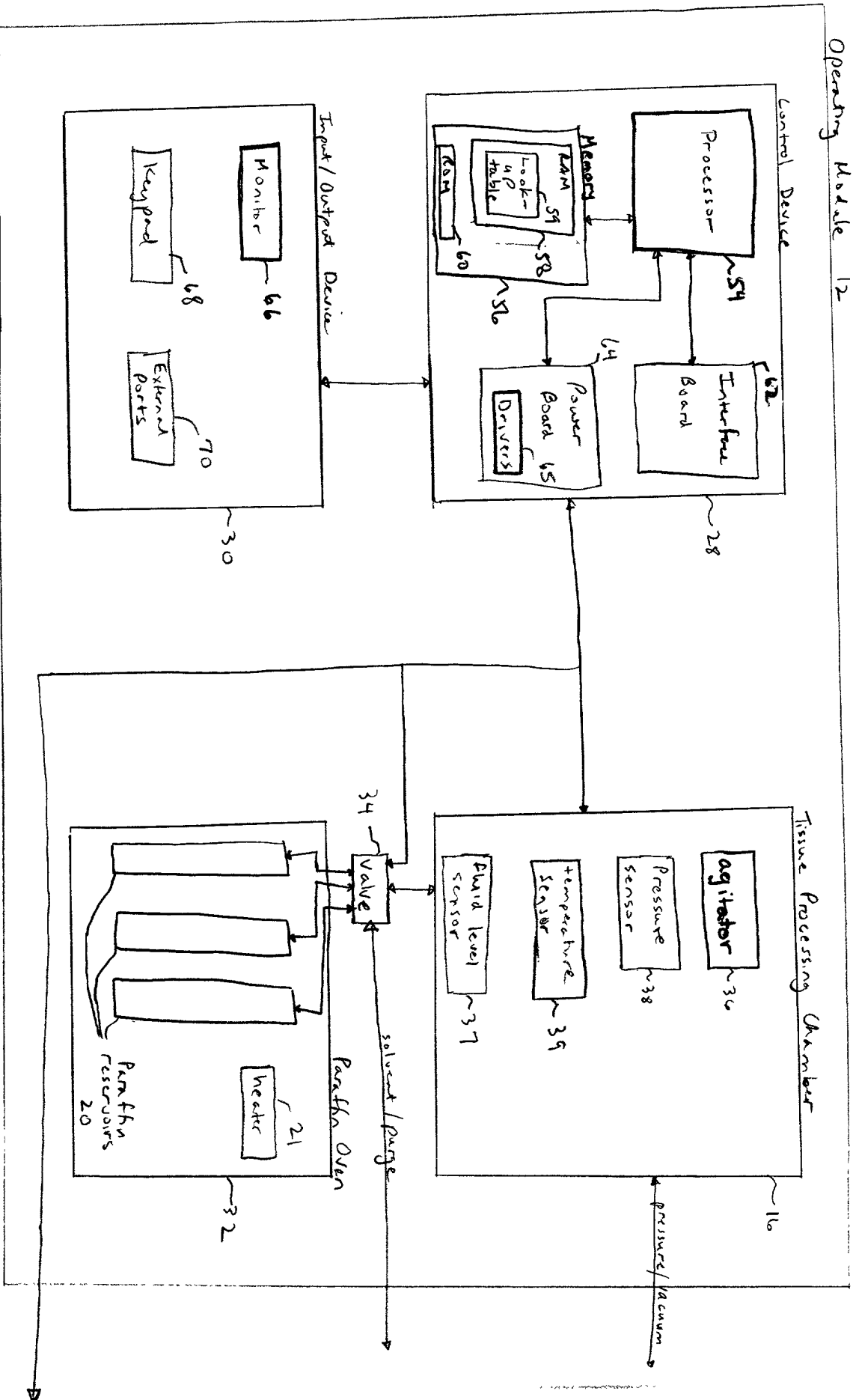


Figure 2a

00040000 00040000

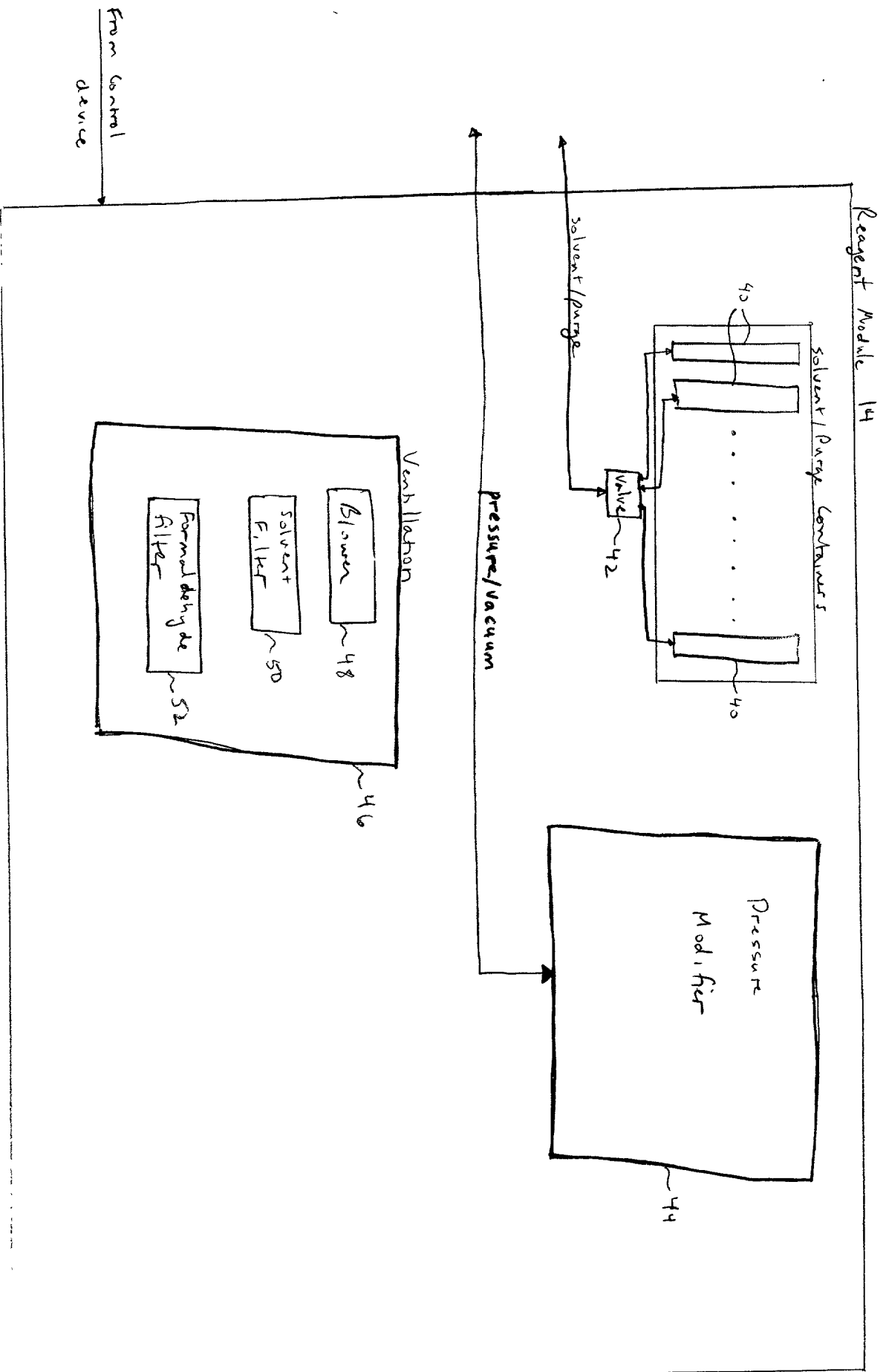


Figure 2b


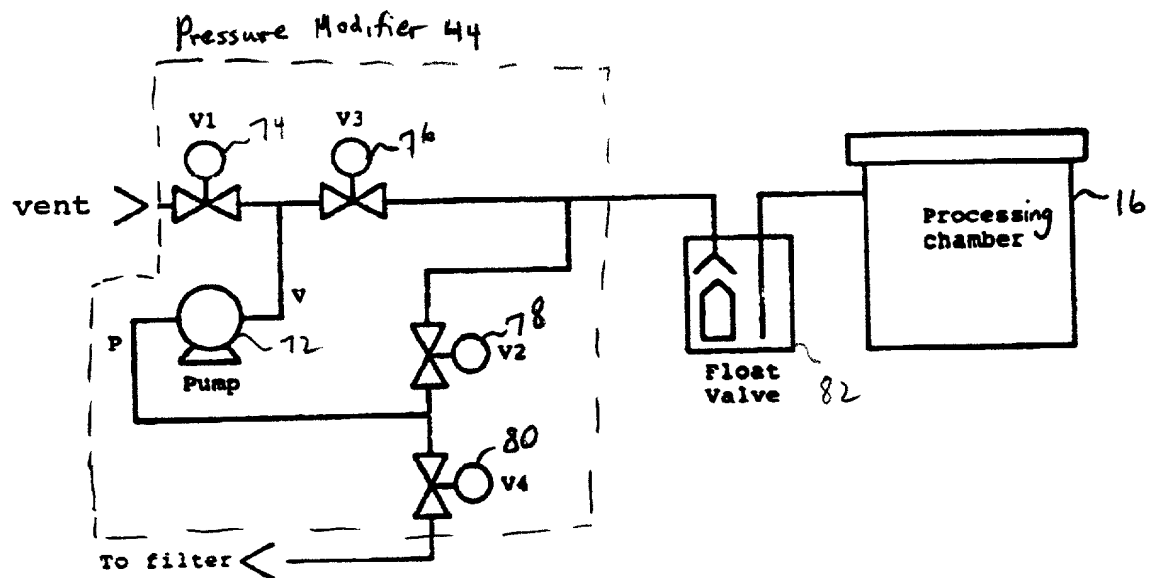


FIGURE 3

2025 RELEASE UNDER E.O. 14176

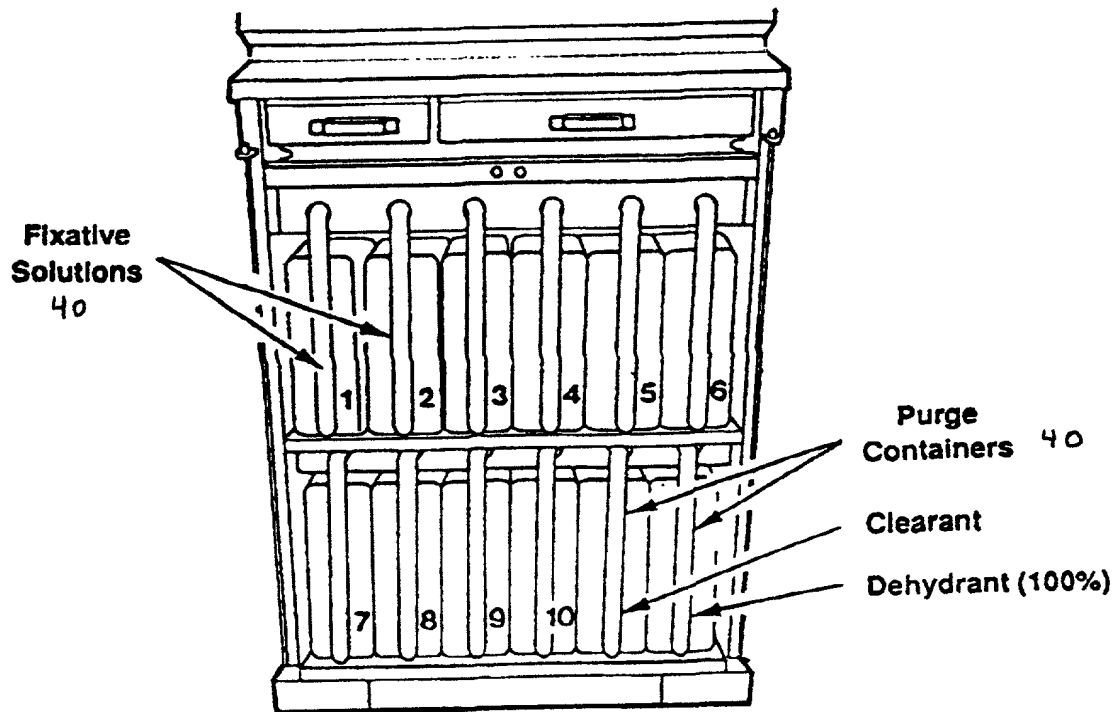

Solvent Container Positions in the Reagent Module

FIGURE 4


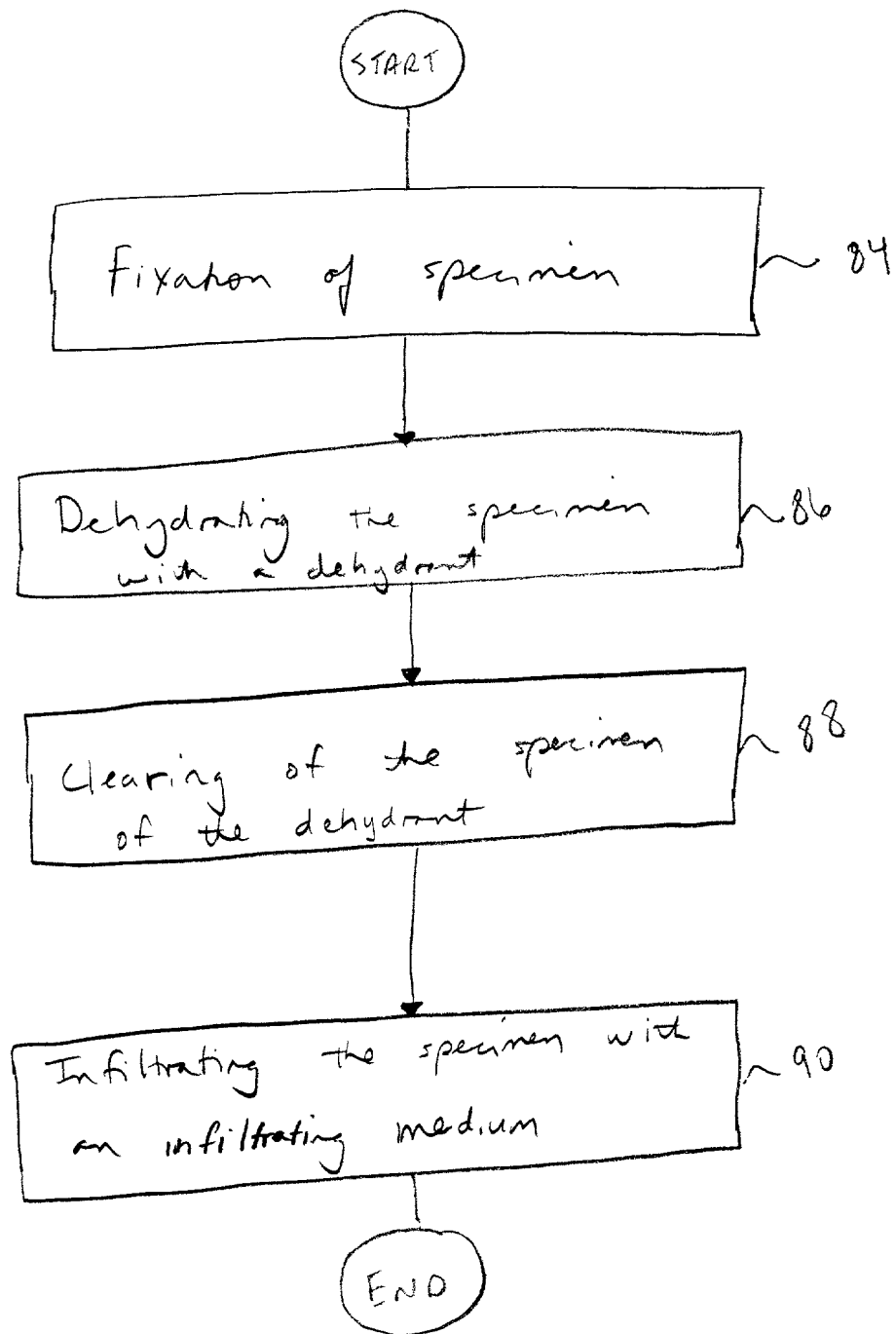


FIGURE 5

FIGURE 6a

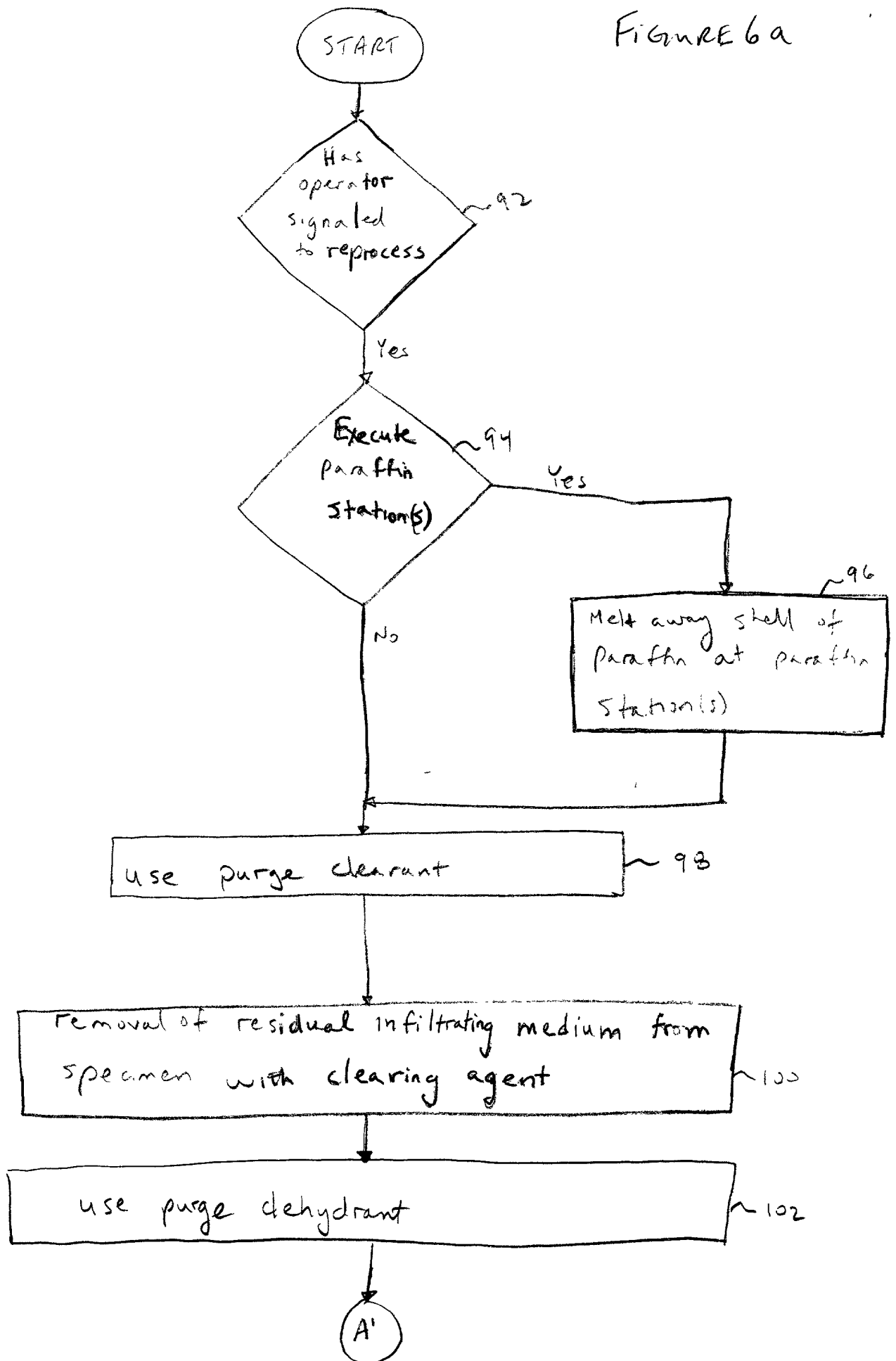

A'

Removal of residual clearing agent by saturating specimen with dehydrating agent 104

Removal of dehydrating agent with aqueous fluid 106

Has operator indicated to process specimen 108

Yes


refixed specimen is dehydrated using dehydrant 110

dehydrant in specimen is replaced using clearing agent 112

clearing agent is replaced using infiltrating material 114

end

FIGURE 6b

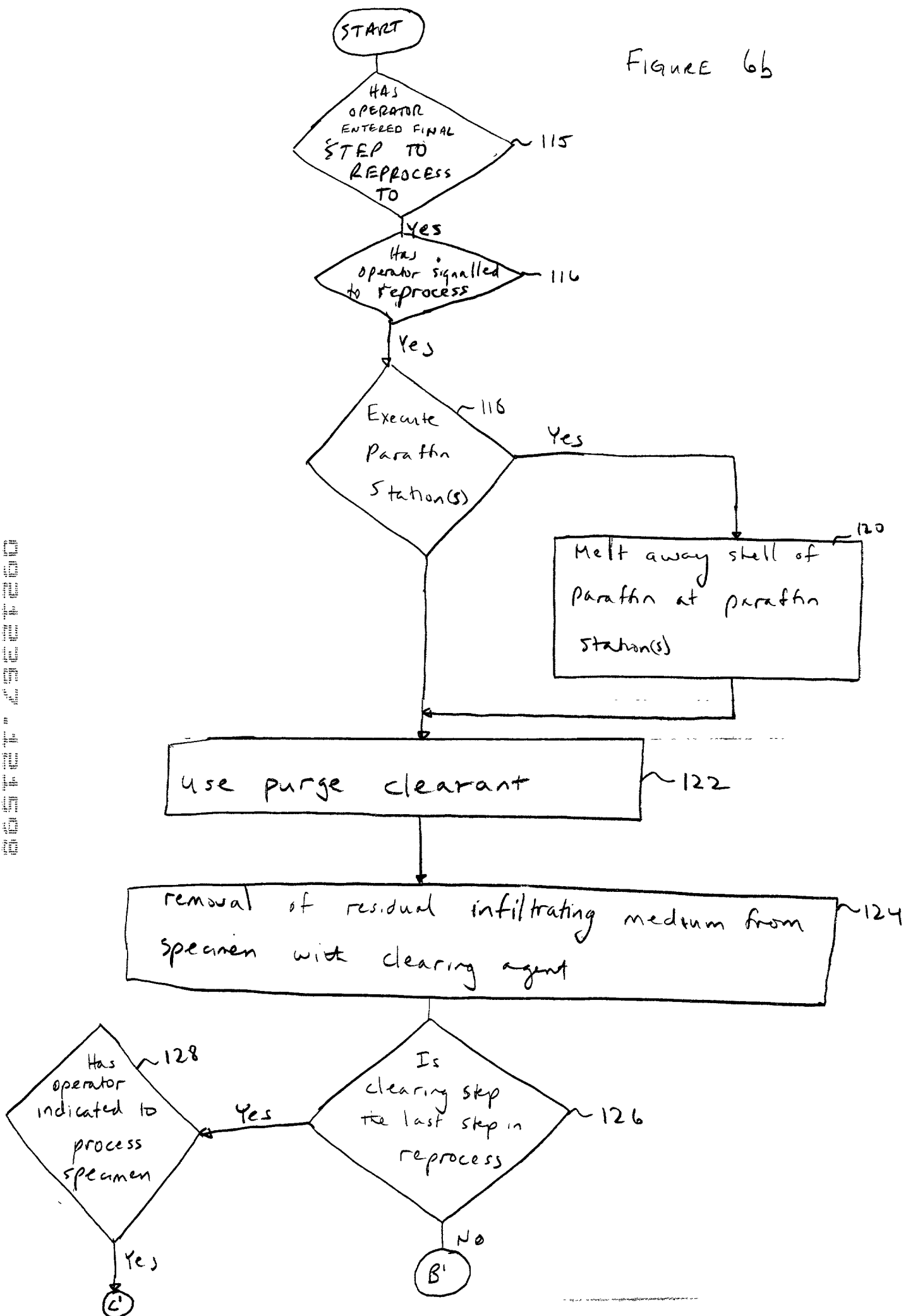

C'

B'

use purge dehydrant ~ 130

removal of residual clearing agent by saturating specimen with dehydrating agent

132

Has operator indicated to process specimen

136

Yes

Yes

Is dehydrating step the last step in reprocess

134

No

Removal of dehydrating agent with aqueous fluid

138

Has operator indicated to process specimen

140

Yes

refixed specimen is dehydrated using dehydrant

142

dehydrant in specimen is replaced with clearing agent

144

clearing agent is replaced using infiltrating material

146

END

Flowchart illustrating the process of specimen reprocessing, including steps for dehydrating, clearing, and infiltrating.